



English Edition 



REACTORS & ACCESSORIES

Version 4 / Edition Summer 2025

Schmizo AG - swiss scientific glass

75 years of excellent glass products for laboratories

Schmizo AG is your technical glassblowing company for custom-made glass products. Founded 75 years ago in Oftringen, we are an owner-managed company with decades of experience and craftsmanship, manufacturing individual glass solutions for a wide range of applications. Whether in industry, science or for custom-made products - we combine precision and creativity with the highest quality.

Above-average quality standards

Our strength lies in our above-average quality standards in the development and manufacture of customer-specific glass apparatus and laboratory reactors. Direct contact with our customers in research, pharmaceuticals and industry is extremely important to us.

Outstanding team of employees

Our foundation consists of a broadly skilled team of 20 employees (10 trained glass apparatus makers, 5 technical staff, and 5 administrative staff). We train our own apprentices, but can also rely on specialists with over 30 years of professional experience.

Flexible partner for borosilicate glass

As a flexible partner, we not only have expertise in the craft of technical glassblowing, but also have an extensive range of machinery for the mechanical precision processing of borosilicate glass.



The Schmizo Reactor Catalogue

This reactor catalogue presents a selection from our wide range of reactors. You will also find the standard accessories required to configure a reactor system to suit your individual needs.

Our technical advisors will be happy to assist you in selecting the right components to ensure an optimally functioning system.



Table of Contents

Flat Flange Reactors and Lids, Filter Systems and Reactor Holders

*Clicking on the desired chapter will take you to the corresponding page.
The blue arrow takes you back to the table of contents.*

Schmizo Reactor Systems	4
Flanges & Tempering Ports	8
Unjacketed Reaction Vessels	9
Round Bottom Vessel, Cylindrical Shape	9
Round Bottom Vessel, Wide Shape	10
Reaction Vessel with Bottom Outlet	11
Flat Bottom Vessel, Cylindrical Shape	12
Flat Bottom Vessel, Wide Shape	13
Round Bottom Flask	14
Round Bottom Flask with Bottom Outlet	15
Jacketed Reaction Vessels	16
Jacketed Reaction Vessel with Round Bottom	16
Jacketed Reaction Vessel with Flat Bottom	18
Jacketed Reaction Vessel with Bottom Outlet	20
Jacketed Conical Reaction Vessel with Bottom Outlet	22
Jacketed Sulphonating Flask with Bottom Outlet	23
Triple-Walled Reaction Vessels	24
Vacuum Jacketed Reaction Vessel with Bottom Outlet	24
Vacuum Jacketed Conical Reaction Vessel with Bottom Outlet	25
Column Extender	26
Column Extender	26
Column Extender Jacketed	27
Jacketed Beaker	28
Jacketed Beaker	28
Flat Flange Lids	29
Flat Flange Lid DN 60	29
Flat Flange Lid DN 100	30
Flat Flange Lid DN 120	31
Flat Flange Lid DN 150	32
Flat Flange Lid DN 200	33
Filter Systems	34
Filter Element and Filter Support	34
Reactor with Glass Filter Plate	36
Seals and Holders for Flat Flanges	37
O-Rings & Flat Seals	37
Reactor Holders	38

Table of Contents

Reactors Accessories & Schlenk Glassware

*Clicking on the desired chapter will take you to the corresponding page.
The blue arrow takes you back to the table of contents.*

Bottom Discharge Valves	40
Bottom Discharge Valve SCHMIZO	40
Equipment for Inlet/Outlet	41
Hose Adapters	41
Metal Hoses	43
Stirrers	44
Glass Stirrers	44
Glass Stirrers with PTFE Blades	46
Stainless Steel Stirrers	49
BOLA Stirrers Stainless Steel/PTFE	50
Stirrer Bearings	51
Stirrer Bearings with Glass Shaft Guide	51
Stirrer Bearings with Synthetic Guide	52
Magnetic Couplings	55
Mechanical Seals	56
Stirrer Couplings & Stirrer Motors	57
Stirrer Couplings	57
CAT Stirrer Motors	58
Reactors Accessories	59
Liquid Dividers	59
Gas Inlet Tubes & Vacuum Pieces	60
Thermometer	61
Temperature Probes	62
Probe Cables & Tubes	64
pH and Redox Electrodes	65
Laboratory Racks	66
Table & Wall Racks	66
Floor Standing Racks	67
Accessories	68
Schlenk Glassware	69
Schlenk Lines and Flasks	69
Appendix	70
Products & Services Overview	70
Technical Information	71
Contacts	72

SCHMIZO Reactor Systems with Table Rack


Features

- Volumes from 0.1 to 5 liters
- Available with single jacket and vacuum jacket
- Covers compatible with all reactor sizes
- Temperature range -80°C to $+200^{\circ}\text{C}$
(up to $+250^{\circ}\text{C}$ with special materials)
- Dry running mechanical seal (optionally with ceramic ball bearing for aggressive media). Vacuum proof 1mbar
- Innovative zero dead space bottom discharge valve

Components of a Reactor System with Table Rack:

Nr.	Article	Find on Page
1	CAT stirrer motors	58
2	Stirrer couplings	57
3	Coolers	5
4	Dropping funnels	5
5	Mechanical seals	51 - 56
6	Temperature probes	62 - 63
7	Guiding tubes	64
8	Reactor lids	29 - 33
9	Reactor holders	38
10	Reactors	9 - 25
*	Tempering ports	41 - 43
11	Stirrers	44 - 50
12	Bottom discharge valves	40
13	Table racks	66
14	Collecting trays	68
15	Temperature measuring device	63
16	Probe cables	64

* not shown on this picture

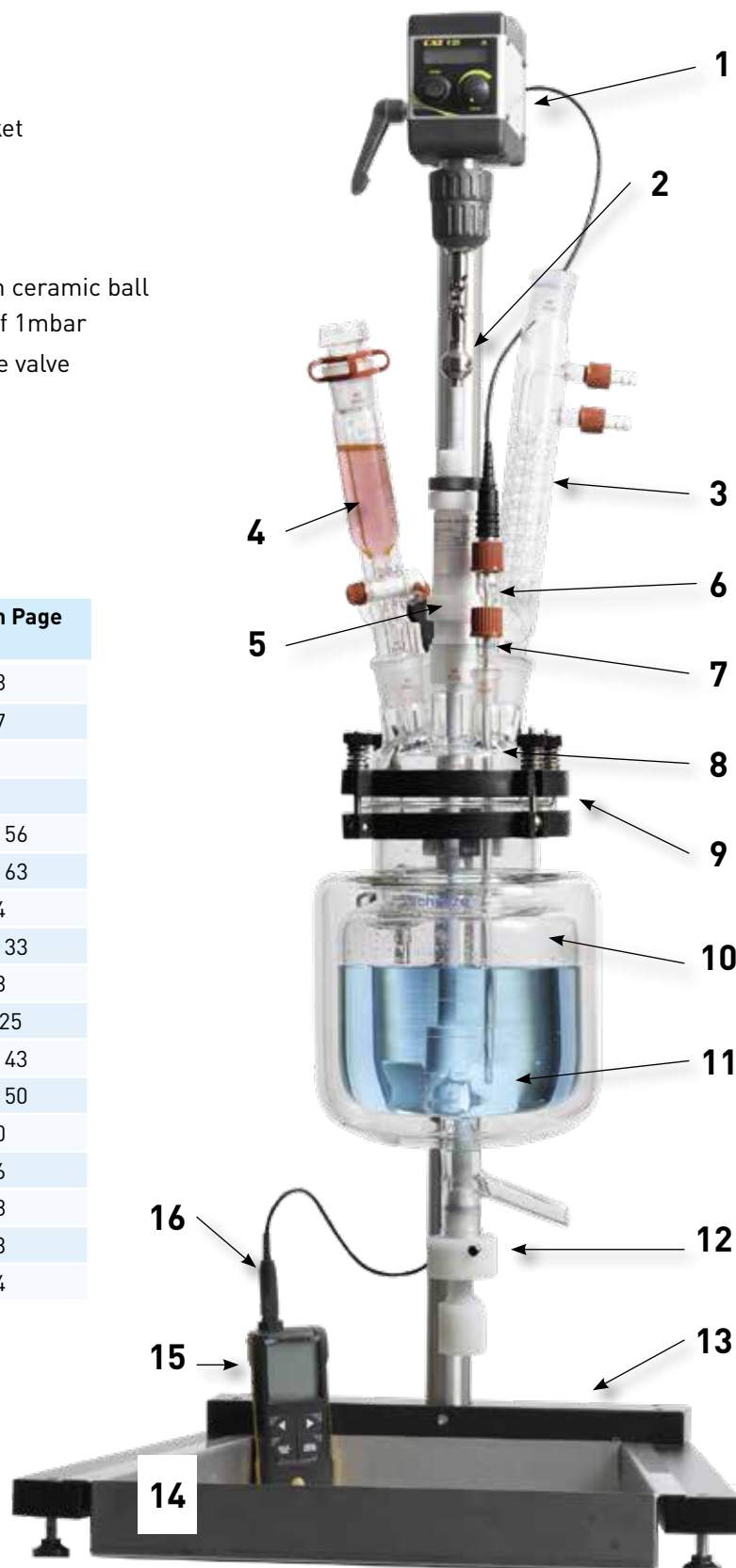

Reactor Accessories

You will find a wide selection of standard coolers, dropping funnels and other ground glass components in the Schmizo **Laboratory Glass and Labware** Catalogue (only as a PDF download on www.schmizo.ch/downloads).

All page references on this page refer to the Laboratory Glass and Labware Catalogue

e.g. Coolers

Article	Page
Ball-Cooler „Allihn“	76
Spiral-Cooler „Graham“	77
„Dimroth“ Cooler	77
„Dimroth“ Cooler with double helix	77
Immersion-Cooler	78
Air-Cooler	78
„Liebig“-Cooler	78
Double area cooler „Davis“	79
Snake-Cooler	80
Ball-Cooler „Allihn“ with NS cone	80
Spiral-Cooler with NS cone	81
Reflux-Cooler „Dimroth“ with NS-Cone	81
Intensiv-Cooler double walled	82
Multiple Spiral-Cooler	82

Photo top right: Dimroth reflux-condenser with NS-cones 29/32 (No. 3 in the illustration on page 4).



e.g. Dropping Funnels without Graduation

Article	Page
Dropping Funnel with glass key	103
Dropping Funnel with PTFE key	104
Dropping Funnel with pressure equilibration, glass key	104
Dropping Funnel with pressure equilibration, PTFE key	105

e.g. Dropping Funnels with Graduation

Article	Page
Dropping Funnel with glass key	105
Dropping Funnel with PTFE key	106
Dropping Funnel with pressure equilibration, glass key	106
Dropping Funnel with pressure equilibration, PTFE key	107

Photo on the right: Graduated dropping funnel with pressure equalization, PTFE key, NSH/NSK 29/32 (No. 4 on the illustration on page 4).



SCHMIZO Reactor-Systems with Floor Standing Rack


Features

- Volumes from 5 to 20 liters
- Available with single jacket and vacuum jacket
- Covers compatible with all reactor sizes
- Temperature range -80°C to $+200^{\circ}\text{C}$
(up to $+250^{\circ}\text{C}$ with special materials)
- Dry running mechanical seal (optionally with ceramic ball bearing for aggressive media). Vacuum proof 1mbar
- Innovative zero dead space bottom discharge valve

Components of a Reactor System with Floor Standing Rack:

Nr.	Artikel	Seiten
1	CAT stirrer motors	58
2	Stirrer couplings	57
3	Coolers	5 / 7
4	Probe cables	64
5	Temperature probes	62 - 63
6	Guiding tubes	64
7	Dropping funnels	5 / 7
8	Mechanical seals	51 - 56
9	Reactor lids	29 - 33
10	Reactor holders	39
11	Reactors	9 - 25
*	Tempering ports	41 - 43
12	Stirrers	44 - 50
13	Bottom discharge valves	40
14	Floor standing racks	67


* nicht shown on this picture

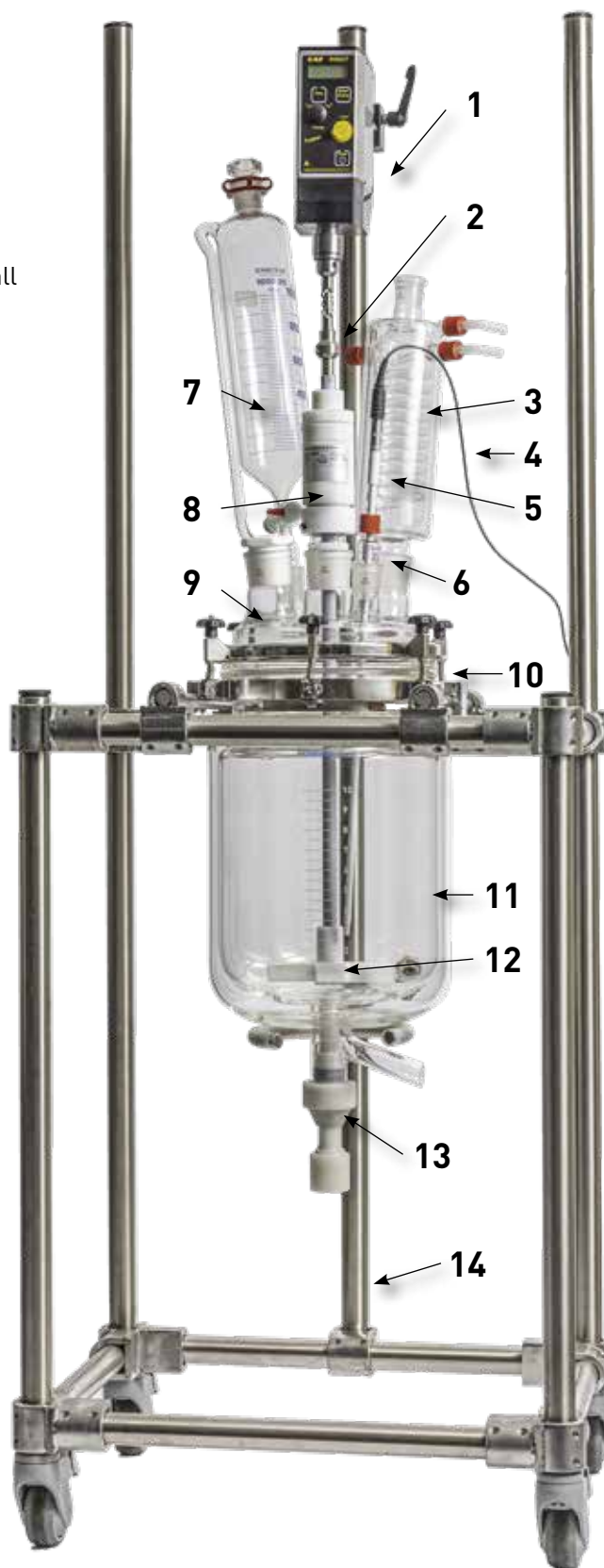

Custom Made Products

In addition to our standard range, we also manufacture laboratory glass components according to your individual requirements.

The example below shows the reactor shown on page 6 with a distillation setup.

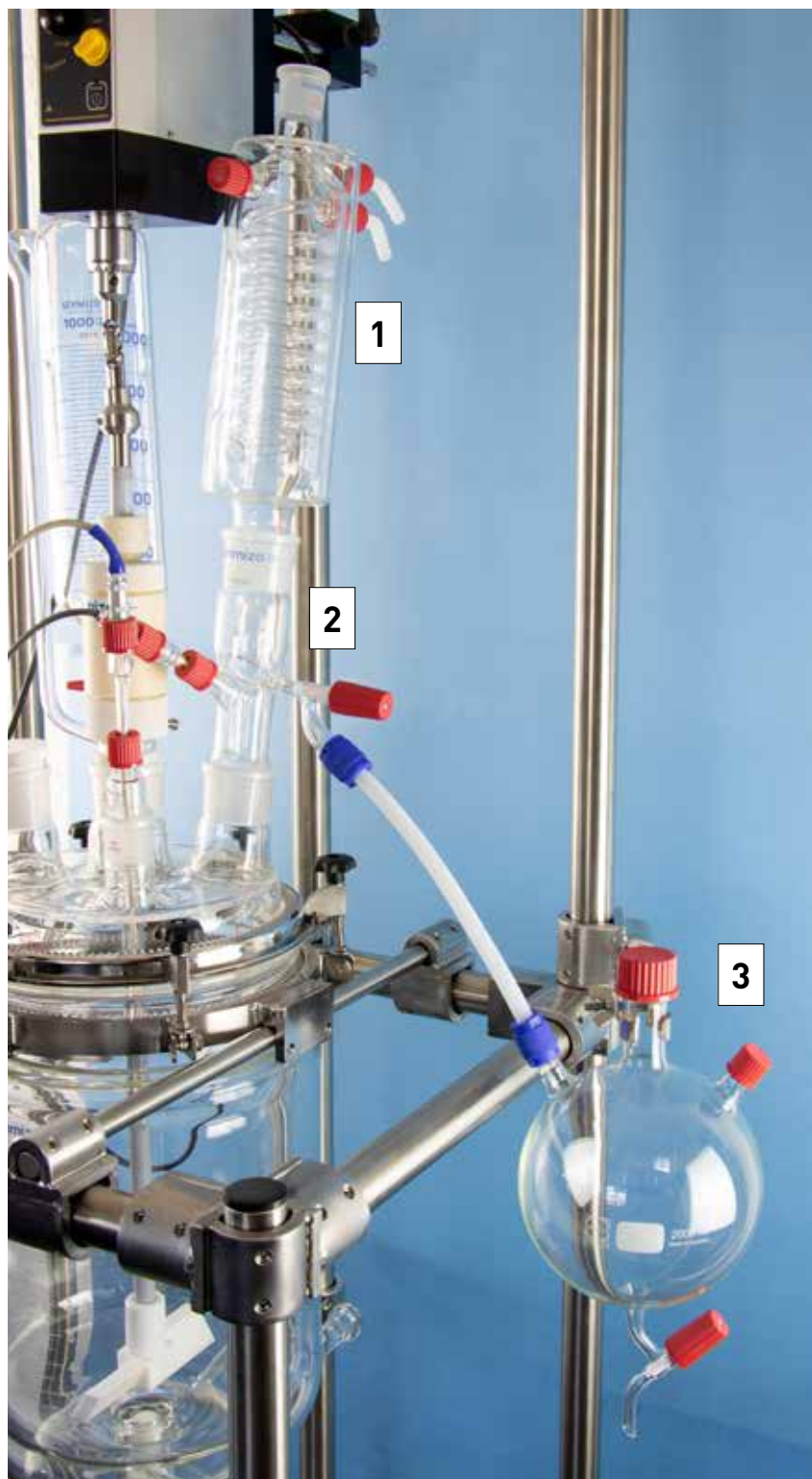

Distillation setup with intensive cooler, liquid divider and receiver



Intensive cooler




Liquid divider




Receiver

Flat Flanges

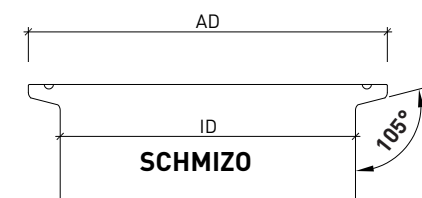
SCHMIZO Flat Flange

DN	ID mm	AD mm
60	68	98
100	105	137
120	120	150
150	150	184

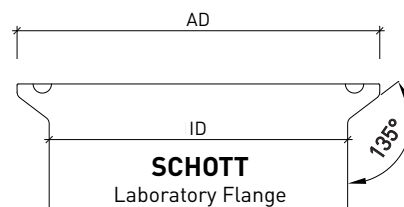
SCHOTT Laboratory Flange


DN	ID mm	AD mm
60	59	100
100	96	138
120	118	158
150	141	184
200	199	242

Tolerances by analogy to DURAN glass tubes.


DN = nominal width , ID = inner diameter, AD = outer diameter

Tempering Ports

AD = outer diameter, ID = inner diameter

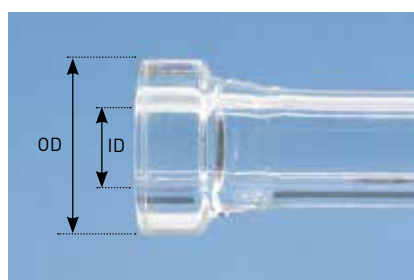

Flat Flange / KF3

DN	ID mm	OD mm
15	15	28.6
25	25	42.2
40	40	57.4
50	50	70

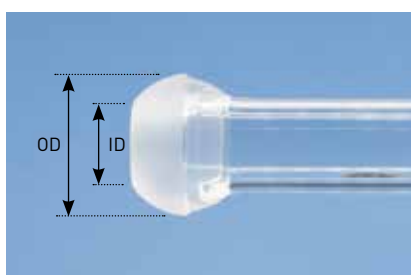
Spherical Joints (KS)

KS	DN mm	OD mm
19/9	9	19
29/15	15	29

Other sizes by request

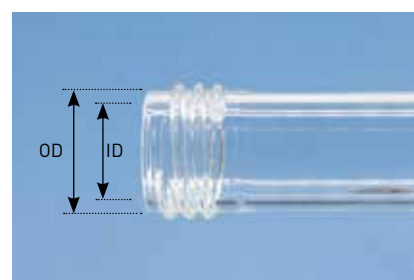

Screwthread Tube

Type	ID mm	OD mm
GL14	9	12
GL18	12.4	16
GL25	18.4	22
SVL15	12.4	16
SVL22	18.4	22

On request also available with hose connections (AD 6 - 12 mm). All other kinds of flanges or connection systems by request.

Round Bottom Vessel, Cylindrical Shape

Reaction vessel, cylindrical shape, round bottom, flat flange SCHMIZO with a groove

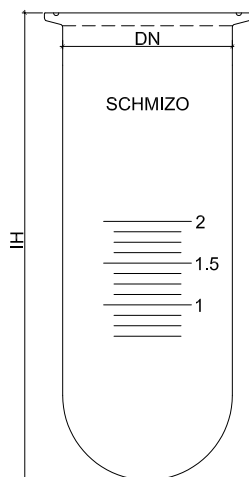

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Cylindrical Shape

The inner diameter is equal to the vessel opening.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
16.100.55	100	60	60	69	-	
16.100.63	200	60	105	69	-	
16.100.65	250	60	120	69	-	
16.100.67	300	60	145	69	-	
16.100.70	350	60	160	69	-	
16.100.73	500	60	200	69	-	
16.110.65	250	100	85	105	-	
16.110.73	500	100	135	105	-	
16.110.79	750	100	175	105	-	
16.110.81	1000	100	205	105	-	
16.110.85	1500	100	255	105	-	
16.110.87	2000	100	315	105	-	
16.120.73	500	120	120	120	-	
16.120.79	750	120	150	120	-	
16.120.81	1000	120	175	120	-	
16.120.85	1500	120	210	120	-	
16.120.87	2000	120	270	120	-	
16.120.88	2500	120	250	120	-	
16.130.81	1000	150	130	150	-	
16.130.85	1500	150	180	150	-	
16.130.87	2000	150	210	150	-	
16.130.88	2500	150	250	150	-	
16.130.89	3000	150	290	150	-	
16.132.94S	5000	200	260	200	-	
16.132.95S	6000	200	305	200	-	
16.132.96S	10000	200	420	200	-	
16.132.97S	15000	200	590	200	-	
16.132.98S	20000	200	770	200	-	

DN 60

DN 100


DN 120

DN 150

DN 200

Round Bottom Vessel, Wide Shape

Reaction vessel, wide shape, round bottom, flat flange SCHMIZO with a groove

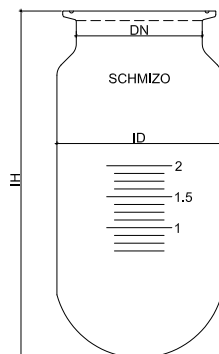

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Wide Shape

The inner diameter of the vessel is bigger than the vessel opening.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
16.600.73	500	60	120	105	-	
16.600.79	750	60	160	105	-	
16.600.81	1000	60	190	105	-	
16.610.85	1500	100	230	120	-	
16.610.87	2000	100	290	120	-	
16.610.88	2500	100	230	160	-	
16.610.89	3000	100	260	160	-	
16.610.92	4000	100	300	160	-	
16.610.94	5000	100	300	180	-	
16.610.95	6000	100	350	180	-	
16.610.96	10000	100	350	240	-	
16.610.97	15000	100	490	240	-	
16.610.98	20000	100	535	260	-	
16.620.89	3000	120	260	160	-	
16.620.92	4000	120	300	160	-	
16.620.94	5000	120	300	180	-	
16.620.95	6000	120	350	180	-	
16.620.96	10000	120	350	240	-	
16.620.97	15000	120	490	240	-	
16.620.98	20000	120	535	260	-	
16.630.89	3000	150	260	160	-	
16.630.92	4000	150	300	160	-	
16.630.94	5000	150	310	180	-	
16.630.95	6000	150	350	180	-	
16.630.96	10000	150	350	240	-	
16.630.97	15000	150	490	240	-	
16.630.98	20000	150	535	260	-	
16.631.94S	5000	200	310	200	-	
16.631.95S	6000	200	345	200	-	
16.631.96S	10000	200	390	240	-	
16.631.97S	15000	200	460	260	-	
16.631.98S	20000	200	550	260	-	

DN 60

DN 100


DN 120

DN 150

DN 200

Reaction Vessel with Bottom Outlet

Reaction vessel, wide shape, round bottom, with bottom outlet opening for discharge valve type SCHMIZO, flat flange SCHMIZO with groove

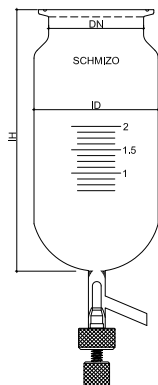

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Discharge Valve SCHMIZO

not included; see page 36.

Side Arm

On request with hose nozzle, spherical joint or GL/SVL thread on the pipe available.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
16.602.74	500	60	120	105	10mm	
16.602.79	750	60	160	105	10mm	
16.602.81	1000	60	190	105	10mm	

DN 60 / 10 mm

16.612.85	1500	100	230	120	10mm	
16.612.87	2000	100	290	120	10mm	
16.612.88	2500	100	230	160	10mm	
16.612.89	3000	100	260	160	10mm	
16.612.92	4000	100	300	160	10mm	
16.612.94	5000	100	300	180	10mm	
16.612.95	6000	100	350	180	10mm	
16.612.96	10000	100	350	240	16mm	
16.612.97	15000	100	490	240	16mm	
16.612.98	20000	100	535	260	16mm	

DN 100 / 10 mm

DN 100 / 16 mm

16.622.89	3000	120	260	160	10mm	
16.622.92	4000	120	300	160	10mm	
16.622.94	5000	120	300	180	10mm	
16.622.95	6000	120	350	180	10mm	
16.622.96	10000	120	350	240	16mm	
16.622.97	15000	120	490	240	16mm	
16.622.98	20000	120	535	260	16mm	

DN 120 / 10 mm

DN 120 / 16 mm

16.632.89	3000	150	230	160	10mm	
16.632.92	4000	150	300	160	10mm	
16.632.94	5000	150	300	180	10mm	
16.632.95	6000	150	350	180	10mm	
16.632.96	10000	150	350	240	16mm	
16.632.97	15000	150	490	240	20mm	
16.632.98	20000	150	535	260	20mm	

DN 150 / 10 mm

DN 150 / 16 mm


DN 150 / 20 mm

16.633.94S	5000	200	310	200	20mm	
16.633.95S	6000	200	345	200	20mm	
16.633.96S	10000	200	390	240	20mm	
16.633.97S	15000	200	460	260	20mm	
16.633.98S	20000	200	550	260	20mm	

DN 200 / 20 mm

Flat Bottom Vessel, Cylindrical Shape

Reaction vessel, cylindrical shape, flat bottom, flat flange SCHMIZO with a groove

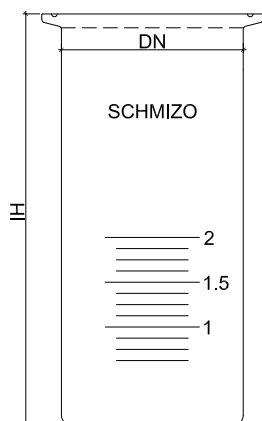

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Cylindrical Shape

The inner diameter is equal to the vessel opening.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
16.200.55	100	60	45	69	-	
16.200.63	200	60	90	69	-	
16.200.65	250	60	105	69	-	
16.200.67	300	60	130	69	-	
16.200.70	350	60	145	69	-	
16.200.73	500	60	195	69	-	
16.200.79	750	60	290	69	-	
16.200.81	1000	60	390	69	-	
16.210.65	250	100	60	105	-	
16.210.73	500	100	110	105	-	
16.210.79	750	100	150	105	-	
16.210.81	1000	100	180	105	-	
16.210.85	1500	100	230	105	-	
16.210.87	2000	100	290	105	-	
16.220.73	500	120	90	120	-	
16.220.79	750	120	120	120	-	
16.220.81	1000	120	145	120	-	
16.220.85	1500	120	190	120	-	
16.220.87	2000	120	240	120	-	
16.220.88	2500	120	250	120	-	
16.230.81	1000	150	90	150	-	
16.230.85	1500	150	140	150	-	
16.230.87	2000	150	170	150	-	
16.230.88	2500	150	210	150	-	
16.230.89	3000	150	250	150	-	
16.232.94S	5000	200	240	200	-	
16.232.95S	6000	200	285	200	-	
16.232.96S	10000	200	400	200	-	
16.232.97S	15000	200	570	200	-	
16.232.98S	20000	200	750	200	-	

DN 60

DN 100


DN 120

DN 150

DN 200

Flat Bottom Vessel, Wide Shape

Reaction vessel, wide shape, flat bottom, flat flange SCHMIZO with a groove

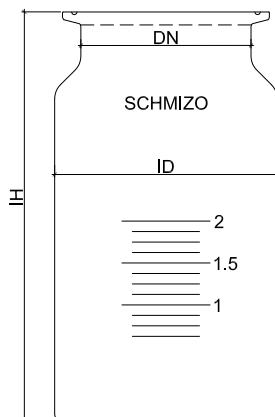

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Wide Shape

The inner diameter of the vessel is bigger than the vessel opening.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
16.700.73	500	60	120	105	-	
16.700.79	750	60	160	105	-	
16.700.81	1000	60	190	105	-	
16.710.85	1500	100	230	120	-	
16.710.87	2000	100	290	120	-	
16.710.88	2500	100	230	160	-	
16.710.89	3000	100	260	160	-	
16.710.92	4000	100	300	160	-	
16.710.94	5000	100	323	180	-	
16.710.96	10000	100	380	220	-	
16.710.97	15000	100	410	260	-	
16.710.98	20000	100	500	280	-	
16.720.89	3000	120	260	160	-	
16.720.92	4000	120	300	160	-	
16.720.94	5000	120	323	180	-	
16.720.96	10000	120	380	220	-	
16.720.97	15000	120	410	260	-	
16.720.98	20000	120	500	280	-	
16.730.89	3000	150	260	160	-	
16.730.92	4000	150	300	160	-	
16.730.94	5000	150	323	180	-	
16.730.96	10000	150	380	220	-	
16.730.97	15000	150	410	260	-	
16.730.98	20000	150	500	280	-	
16.732.96S	10000	200	380	220	-	
16.732.97S	15000	200	410	260	-	
16.732.98S	20000	200	500	280	-	

DN 60

DN 100


DN 120

DN 150

DN 200

Round Bottom Flask

Reaction vessel, round bottom flask, flat flange SCHMIZO with a groove

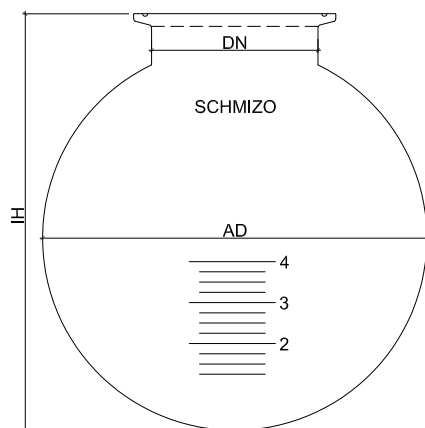

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
16.300.73	500	60	150	105	-	
16.300.81	1000	60	180	131	-	
16.300.87	2000	60	210	166	-	
16.310.81	1000	100	145	131	-	
16.310.87	2000	100	220	166	-	
16.310.89	3000	100	240	185	-	
16.310.92	4000	100	260	206	-	
16.310.94	5000	100	275	223	-	
16.310.95	6000	100	290	236	-	
16.310.96	10000	100	315	279	-	
16.320.87	2000	120	220	166	-	
16.320.89	3000	120	240	185	-	
16.320.92	4000	120	260	206	-	
16.320.94	5000	120	275	223	-	
16.320.95	6000	120	290	236	-	
16.320.96	10000	120	315	279	-	
16.330.89	3000	150	240	185	-	
16.330.92	4000	150	260	206	-	
16.330.94	5000	150	275	223	-	
16.330.95	6000	150	290	236	-	
16.330.96	10000	150	315	279	-	
16.340.96S	10000	200	360	279	-	
16.340.98S	20000	200	420	345	-	

DN 60

DN 100


DN 120

DN 150

DN 200

Round Bottom Flask with Bottom Outlet

Reaction vessel, round bottom flask, with bottom outlet opening for discharge valve type SCHMIZO, flat flange SCHMIZO with groove

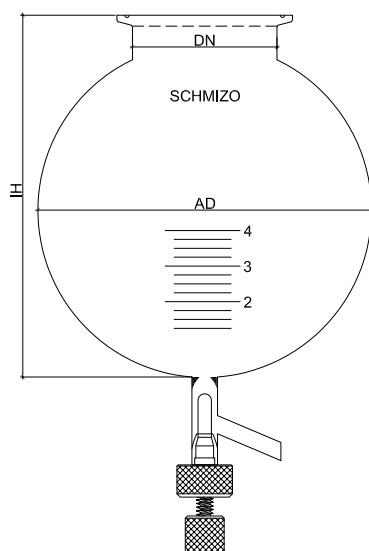

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Discharge-Valve SCHMIZO

not included; see page 36.

Side Arm

On request with hose nozzle, spherical joint or GL/SVL thread on the pipe available.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
16.302.74	500	60	150	105	10mm	
16.302.81	1000	60	180	131	10mm	
16.302.87	2000	60	210	166	10mm	

DN 60 / 10 mm

16.312.81	1000	100	190	131	10mm	
16.312.87	2000	100	220	166	10mm	
16.312.89	3000	100	240	185	10mm	
16.312.92	4000	100	260	206	10mm	
16.312.94	5000	100	275	223	10mm	
16.312.95	6000	100	290	236	10mm	
16.312.96	10000	100	340	279	16mm	

DN 100 / 10 mm

DN 100 / 16 mm

16.322.87	2000	120	220	166	10mm	
16.322.89	3000	120	240	185	10mm	
16.322.92	4000	120	260	206	10mm	
16.322.94	5000	120	275	223	10mm	
16.322.95	6000	120	290	236	10mm	
16.322.96	10000	120	340	279	16mm	

DN 120 / 16 mm

DN 120 / 16 mm

16.332.89	3000	150	240	185	10mm	
16.332.92	4000	150	260	206	10mm	
16.332.94	5000	150	275	223	10mm	
16.332.95	6000	150	290	236	10mm	
16.332.96	10000	150	340	279	16mm	

DN 150 / 10 mm


DN 150 / 16 mm

16.342.96S	10000	200	360	279	20mm	
16.342.98S	20000	200	420	345	20mm	

DN 200 / 20 mm

Jacketed Reaction Vessel with Round Bottom

Reaction vessel, jacketed, flat flange SCHMIZO with a groove, round bottom, inlet/outlet flat flange DN 15

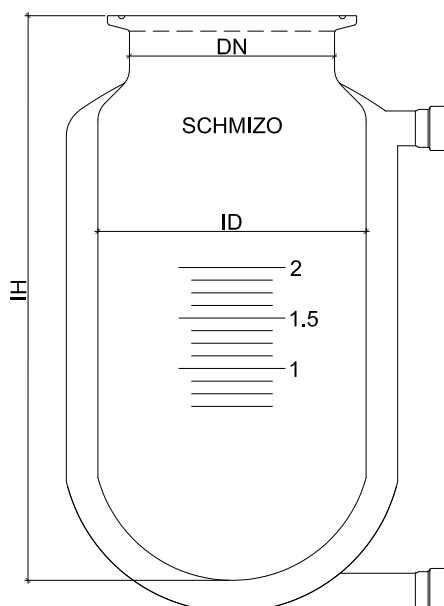

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Inlet/Outlet

Flat flange DN 15

On request also available with spherical joints, GL or SVL threads.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
17.300.55	100	60	95	65	-	
17.300.65	250	60	155	65	-	
17.300.70	350	60	195	65	-	
17.300.73	500	60	250	65	-	
17.300.74	500	60	190	85	-	
17.300.79	750	60	185	105	-	
17.300.81	1000	60	210	105	-	
17.310.73	500	100	150	105	-	
17.310.79	750	100	190	105	-	
17.310.81	1000	100	210	105	-	
17.310.85	1500	100	230	120	-	
17.310.87	2000	100	205	160	-	
17.310.88	2500	100	230	160	-	
17.310.89	3000	100	260	160	-	
17.310.92	4000	100	300	160	-	
17.310.94	5000	100	300	180	-	
17.310.95	6000	100	350	180	-	
17.310.96	10000	100	350	240	-	
17.310.97	15000	100	490	240	-	
17.310.98	20000	100	535	260	-	

DN 60

DN 100


Jacketed Reaction Vessel with Round Bottom (Continuation)

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
17.320.81	1000	120	180	120	-	
17.320.85	1500	120	230	120	-	
17.320.87	2000	120	290	120	-	
17.320.88	2500	120	230	160	-	
17.320.89	3000	120	260	160	-	
17.320.92	4000	120	300	160	-	
17.320.94	5000	120	300	180	-	
17.320.95	6000	120	350	180	-	
17.320.96	10000	120	350	240	-	
17.320.97	15000	120	490	240	-	
17.320.98	20000	120	535	260	-	
17.330.88	2500	150	230	160	-	
17.330.89	3000	150	260	160	-	
17.330.92	4000	150	300	160	-	
17.330.94	5000	150	310	180	-	
17.330.95	6000	150	350	180	-	
17.330.96	10000	150	350	240	-	
17.330.97	15000	150	490	240	-	
17.330.98	20000	150	535	260	-	
17.331.94S	5000	200	310	200	-	
17.331.95S	6000	200	345	200	-	
17.331.96S	10000	200	470	200	-	
17.332.96S	10000	200	390	240	-	
17.332.97S	15000	200	460	260	-	
17.332.98S	20000	200	550	260	-	

DN 120


DN 150

DN 200

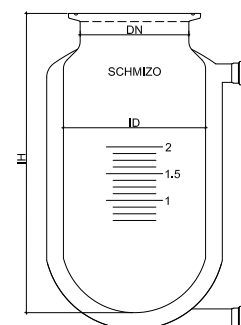

Jacketed Reaction Vessel with Flat Bottom

Reaction vessel, jacketed, flat flange SCHMIZO with a groove, flat bottom, inlet/outlet flat flange DN 15

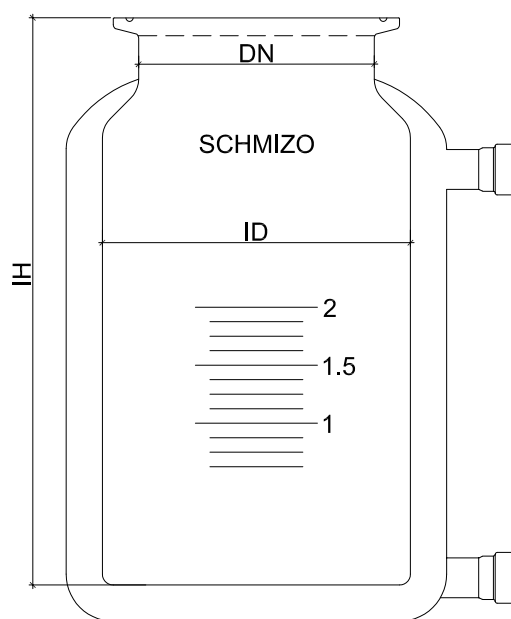

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Inlet/Outlet

Flat flange DN 15

On request also available with spherical joints, GL or SVL threads.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
17.200.55	100	60	85	65	-	
17.200.65	250	60	145	65	-	
17.200.70	350	60	185	65	-	
17.200.73	500	60	240	65	-	
17.200.74	500	60	185	85	-	
17.200.79	750	60	165	105	-	
17.200.81	1000	60	190	105	-	
17.210.73	500	100	125	105	-	
17.210.79	750	100	165	105	-	
17.210.81	1000	100	200	105	-	
17.210.85	1500	100	230	120	-	
17.210.87	2000	100	175	160	-	
17.210.88	2500	100	205	160	-	
17.210.89	3000	100	230	160	-	
17.210.92	4000	100	290	160	-	
17.210.94	5000	100	285	180	-	
17.210.95	6000	100	335	200	-	
17.210.96	10000	100	380	220	-	
17.210.97	15000	100	410	260	-	
17.210.98	20000	100	500	280	-	

DN 60

DN 100

Jacketed Reaction Vessel with Flat Bottom (Continuation)

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
17.220.81	1000	120	170	120	-	
17.220.85	1500	120	220	120	-	
17.220.87	2000	120	275	120	-	
17.220.88	2500	120	205	160	-	
17.220.89	3000	120	230	160	-	
17.220.92	4000	120	290	160	-	
17.220.94	5000	120	285	180	-	
17.220.95	6000	120	335	200	-	
17.220.96	10000	120	380	220	-	
17.220.97	15000	120	410	260	-	
17.220.98	20000	120	500	280	-	
17.230.88	2500	150	200	160	-	
17.230.89	3000	150	230	160	-	
17.230.92	4000	150	290	160	-	
17.230.94	5000	150	290	180	-	
17.230.95	6000	150	335	200	-	
17.230.96	10000	150	380	220	-	
17.230.97	15000	150	410	260	-	
17.230.98	20000	150	500	280	-	
17.231.94S	5000	200	260	180	-	
17.231.95S	6000	200	335	200	-	
17.231.96S	10000	200	450	220	-	
17.232.96S	10000	200	380	240	-	
17.232.97S	15000	200	410	260	-	
17.232.98S	20000	200	500	280	-	


DN 120

DN 150

DN 200

Jacketed Reaction Vessel with Bottom Outlet

Reaction vessel, jacketed, round bottom, with bottom outlet opening for discharge valve type SCHMIZO, flat flange SCHMIZO with groove, inlet/outlet flat flange DN 15

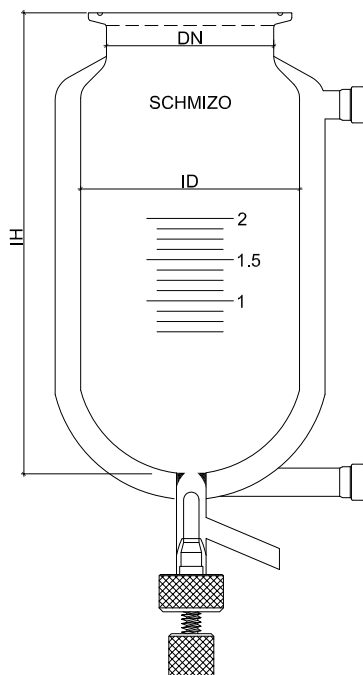

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with **inlet/outlet** flange.

Flat flange DN 15

On request also available with spherical joints, GL or SVL threads

Discharge-Valve SCHMIZO

not included; see page 36.

Side Arm

On request with hose nozzle, spherical joint or GL/SVL thread on the pipe available.

Article No.	Volume ml	DN	IH mm	ID mm	Valve bore	
17.105.55	100	60	95	65	10mm	
17.105.65	250	60	155	65	10mm	
17.105.70	350	60	195	65	10mm	
17.105.73	500	60	250	65	10mm	
17.105.74	500	60	190	85	10mm	
17.105.79	750	60	185	105	10mm	
17.105.81	1'000	60	210	105	10mm	

DN 60 / 10 mm

17.110.73	500	100	150	105	10mm	
17.110.79	750	100	190	105	10mm	
17.110.81	1'000	100	210	105	10mm	
17.110.85	1'500	100	210	120	10mm	
17.110.86	2'000	100	280	120	10mm	
17.110.87	2'000	100	205	160	10mm	
17.110.88	2'500	100	230	160	10mm	
17.110.89	3'000	100	260	160	10mm	
17.110.92	4'000	100	290	160	10mm	
17.110.94	5'000	100	290	180	10mm	
17.110.95	6'000	100	335	180	10mm	
17.110.96	10'000	100	340	240	16mm	
17.110.97	15'000	100	450	240	16mm	
17.110.98	20'000	100	520	260	16mm	

DN 100 / 10 mm

DN 100 / 16 mm

Jacketed Reaction Vessel with Bottom Outlet (Continuation)

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
17.120.81	1'000	120	180	120	10mm	
17.120.85	1'500	120	220	120	10mm	
17.120.87	2'000	120	280	120	10mm	
17.120.88	2'500	120	230	160	10mm	
17.120.89	3'000	120	260	160	10mm	
17.120.92	4'000	120	290	160	10mm	
17.120.94	5'000	120	310	180	10mm	
17.120.95	6'000	120	335	180	10mm	
17.120.96	10'000	120	350	240	16mm	
17.120.97	15'000	120	450	240	16mm	
17.120.98	20'000	120	520	260	16mm	
17.130.88	2'500	150	230	160	10mm	
17.130.89	3'000	150	260	160	10mm	
17.130.92	4'000	150	290	160	10mm	
17.130.94	5'000	150	310	180	10mm	
17.130.95	6'000	150	335	180	10mm	
17.130.96	10'000	150	350	240	16mm	
17.130.97	15'000	150	450	240	16mm	
17.130.98	20'000	150	520	260	16mm	
17.132.94S	5'000	200	260	200	20mm	
17.132.95S	6'000	200	360	200	20mm	
17.132.96S	10'000	200	480	200	20mm	
17.133.96S	10'000	200	360	240	20mm	
17.133.97S	15'000	200	460	260	20mm	
17.133.98S	20'000	200	550	260	20mm	

DN 120 / 10 mm


DN 120 / 16 mm

DN 150 / 10 mm

DN 150 / 16 mm

DN 200 / 20 mm


On request all articles are also available with valve bore 20 mm.

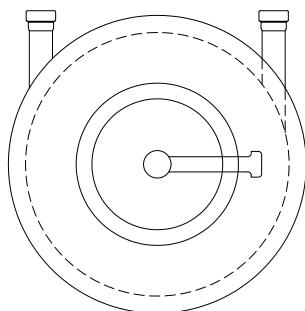

Top view: Inlet/outlet backwards, side arm rightwards. On request we can customize this arrangement.

Jacketed Reaction Vessel, Conical Shape, with Bottom Outlet

Reaction vessel, jacketed, conical shape, with bottom outlet opening for discharge valve type SCHMIZO, flat flange SCHMIZO with groove, inlet/outlet flat flange DN 15

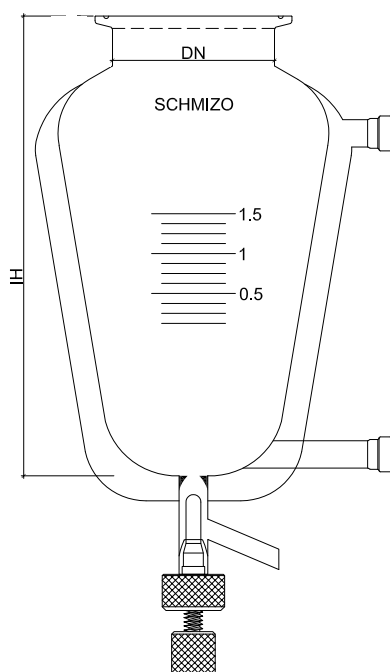

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Inlet/Outlet

Flat flange DN 15

On request also available with spherical joints, GL or SVL threads

Discharge-Valve SCHMIZO

not included; see page 36.

Side Arm


On request with hose nozzle, spherical joint or GL/SVL thread on the pipe available.

Article No.	Volume ml	DN	IH mm	ID mm	Valve bore	
17.500.73	500	100	285		10mm	
17.500.81	1000	100	290		10mm	
17.501.81	1000	120	340		10mm	
17.501.87	2000	120	340		10mm	
17.502.87	2000	150	375		10mm	
17.502.89	3000	150	520		10mm	

DN 100 / 10 mm

DN 120 / 10 mm


DN 150 / 10 mm

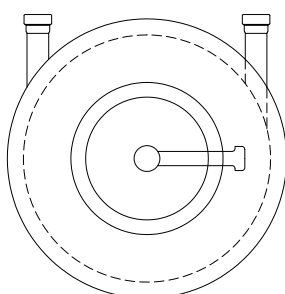

Top view: Inlet/outlet backwards, side arm rightwards. On request we can customize this arrangement.

Jacketed Sulphonating Flask with Bottom Outlet

Jacketed sulphonating flask, with bottom outlet opening for discharge valve type SCHMIZO, inlet/outlet flat flange DN 15

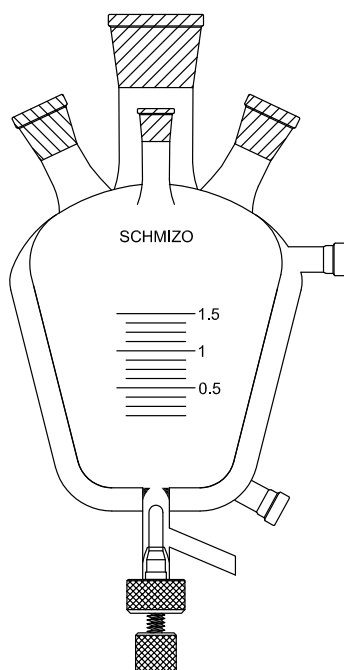

Inlet/Outlet

Flat flange DN 15

On request also available with spherical joints, GL or SVL threads.



Discharge-Valve SCHMIZO

not included; see page 36.

Side Arm

On request with hose nozzle, spherical joint or GL/SVL thread on the pipe available.


Article No.	Volume ml	Central NS	Lateral NS	Lateral NS	Lateral NS	
17.510.72	350	45	29	14	-	
17.510.73	500	45	29	29	-	
17.510.79	750	45	29	29	-	
17.510.85	1500	45	29	29	14	
17.510.89	2500	45	29	29	14	

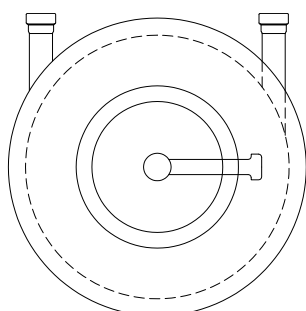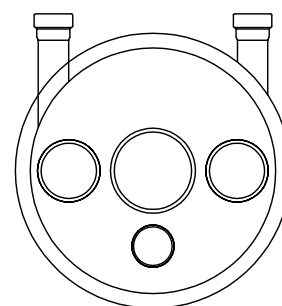

Top view: Inlet/outlet backwards, side arm rightwards. On request we can customize this arrangement.

Vacuum Jacketed Reaction Vessel with Bottom Outlet

Vacuum-triple-walled reaction vessel, round bottom, with bottom outlet opening for discharge valve type SCHMIZO, flat flange SCHMIZO with groove, inlet/outlet flat flange DN 15

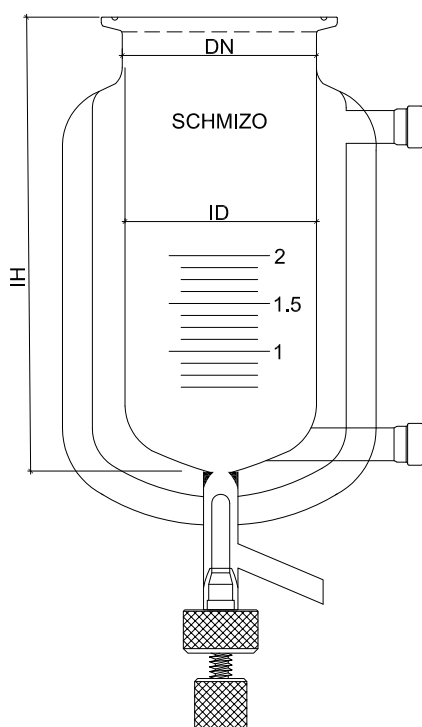

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Inlet/Outlet

Flat flange DN 15

On request also available with spherical joints, GL or SVL threads

Discharge-Valve SCHMIZO

not included; see page 36.

Side Arm

On request with hose nozzle, spherical joint or GL/SVL thread on the pipe available.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
17.520.73	500	60	250	65	10mm	
17.521.73	500	100	190	105	10mm	
17.521.81	1000	100	250	105	10mm	
17.521.87	2000	100	265	150	10mm	
17.521.89	3000	100	320	150	10mm	
17.521.92	4000	100	380	150	10mm	
17.522.81	1000	120	180	120	10mm	
17.522.87	2000	120	290	120	10mm	
17.522.89	3000	120	330	150	10mm	
17.522.92	4000	120	390	150	10mm	
17.522.94	5000	120	360	180	10mm	
17.522.95	6000	120	400	180	10mm	
17.523.87	2000	150	240	150	10mm	
17.523.89	3000	150	285	150	10mm	
17.523.92	4000	150	330	150	10mm	
17.523.94	5000	150	350	180	10mm	
17.523.95	6000	150	370	180	10mm	

— DN 60 / 10 mm


DN 100 / 10 mm

DN 120 / 10 mm

DN 150 / 10 mm

Vacuum Jacketed Reaction Vessel, Conical Shape, with Bottom Outlet

Vacuum-triple-walled reaction vessel, conical shape, with bottom outlet opening for discharge valve type SCHMIZO, flat flange SCHMIZO with groove, inlet/outlet flat flange DN 15

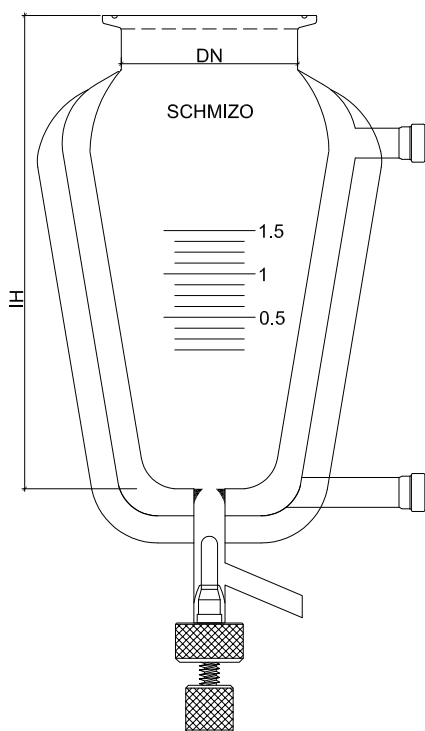

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Inlet/Outlet

Flat flange DN 15

On request also available with spherical joints, GL or SVL threads

Discharge-Valve SCHMIZO

not included; see page 36.


Side Arm

On request with hose nozzle, spherical joint or GL/SVL thread on the pipe available.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
17.530.73	500	100	285		10mm	
17.530.81	1000	100	290		10mm	
17.531.81	1000	100	340		10mm	
17.531.87	2000	100	340		10mm	
17.532.87	2000	120	375		10mm	
17.532.89	3000	120	520		10mm	

DN 100 / 10 mm


DN 120 / 10 mm

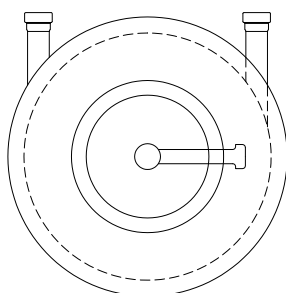

Top view: Inlet/outlet backwards, side arm rightwards. On request we can customize this arrangement.

Column Extender

with two flat flanges type SCHMIZO: 1x with groove, 1x without groove

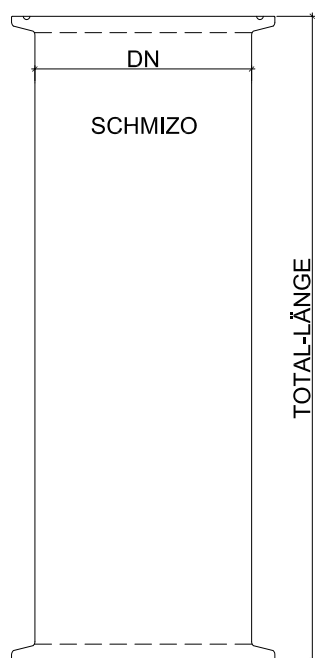

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Article No.	Length mm	DN	
17.540.55	250	60	
17.540.65	400	60	
17.542.63	300	100	
17.542.67	500	100	
17.544.63	300	120	
17.544.65	400	120	
17.544.67	500	120	
17.546.65	400	150	
17.546.67	500	150	
17.548.65S	400	200	
17.548.67S	500	200	
17.548.68S	600	200	
17.548.69S	750	200	

DN 60

DN 100


DN 120

DN 150

DN 200

Column Extender Jacketed

with two flat flanges type SCHMIZO: 1x with groove, 1x without groove;
inlet/outlet flat flange DN 15

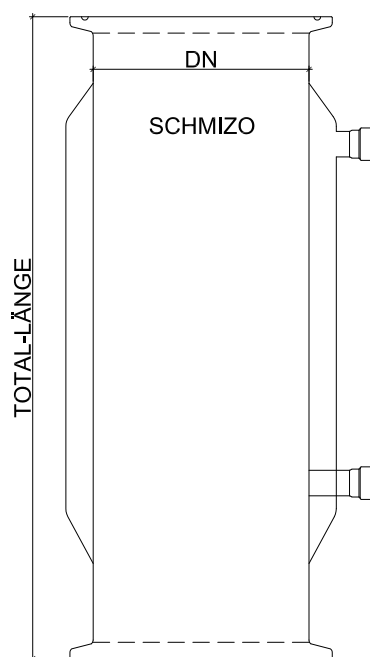

SCHMIZO Flat Flange

DN = Diameter nominal

All reaction vessels are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number. Size DN 200 is only available with SCHOTT-flange.

Inlet/Outlet

Flat flange DN 15

On request also available with spherical joints, GL or SVL threads.

Article No.	Length mm	DN	
17.550.55	250	60	
17.550.65	400	60	
17.552.63	300	100	
17.552.67	500	100	
17.554.63	300	120	
17.554.65	400	120	
17.554.67	500	120	
17.556.65	400	150	
17.556.67	500	150	
17.558.65S	400	200	
17.558.67S	500	200	
17.558.68S	600	200	

DN 60

DN 100


DN 120

DN 150

DN 200

Jacketed Beaker

Jacketed beaker with flat bottom, inlet/outlet flat flange DN 15

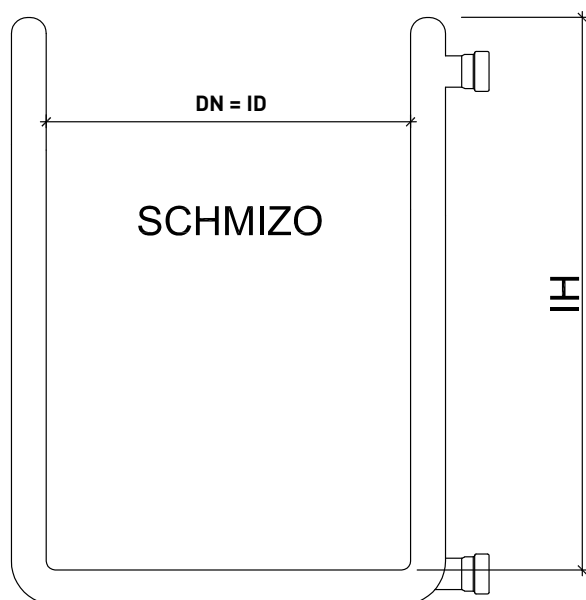

Inlet/Outlet

Flat flange DN 15

On request also available with spherical joints, GL or SVL threads.

Article No.	Volume ml	DN	IH mm	ID mm	Valve Bore	
17.515.55	100	60	60	60	-	
17.515.65	250	60	120	60	-	
17.515.73	500	60	210	60	-	
17.516.81	1000	100	170	100	-	
17.516.87	2000	100	290	100	-	
17.517.87	2000	120	205	120	-	
17.517.88	3000	120	295	120	-	
17.518.89	3000	150	205	150	-	
17.518.92	4000	150	260	150	-	
17.518.94	5000	150	365	150	-	
17.518.95	6000	150	390	150	-	

DN / ID 60


DN / ID 100

DN / ID 120

DN / ID 150

Flat Flange Lid DN 60

Flat Flange SCHMIZO without a groove

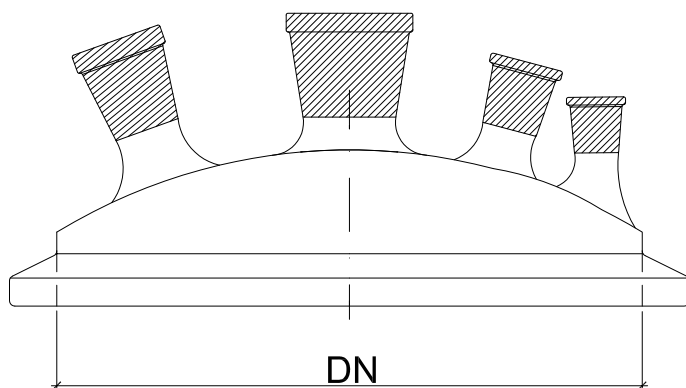

SCHMIZO Flat Flange

DN = Diameter nominal

All lids are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number.

Article No.	DN	Central NS	Lateral Vertical NS	Lateral Angled NS	Lateral Angled NS	Lateral Angled NS	
17.500.00	60						
17.501.10	60	14					
17.501.25	60	29					
17.501.35	60	45					
17.501.45	60	71					
17.502.10	60	14	14				
17.502.25	60	29	14				
17.502.30	60	29	29				
17.503.20	60	29	14	14			
17.503.25	60	29	14	29			
17.504.20	60	29	14	14	14		
17.504.25	60	29	14	14	29		
17.504.30	60	29	14	29	29		
17.505.50	60	29	14	14	29	14	

Blank Lid

1-Neck

2-Neck

3-Neck


4-Neck

5-Neck

Special designs on request

Flat Flange Lid DN 100

Flat Flange SCHMIZO without a groove

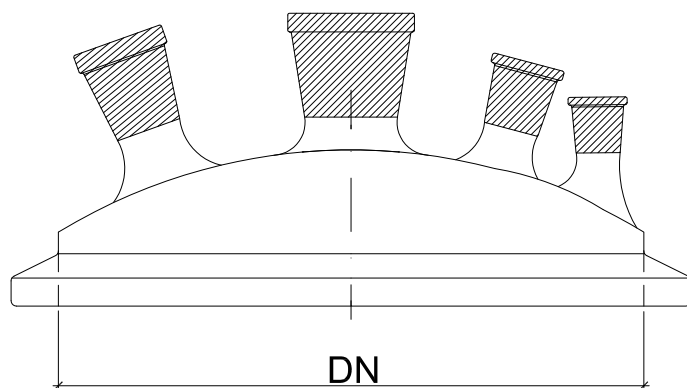

SCHMIZO Flat Flange

DN = Diameter nominal

All lids are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number.

Article No.	DN	Central NS	Lateral Vertical NS	Lateral Angled NS	Lateral Angled NS	Lateral Angled NS	
17.600.00	100						
17.601.10	100	14					
17.601.25	100	29					
17.601.35	100	45					
17.601.45	100	71					
17.602.10	100	14	14				
17.602.25	100	29	14				
17.602.30	100	29	29				
17.602.35	100	45	14				
17.602.40	100	45	29				
17.603.20	100	29	14	14			
17.603.25	100	29	14	29			
17.603.30	100	45	14	14			
17.603.40	100	45	14	29			
17.604.20	100	29	14	14	14		
17.604.25	100	29	14	14	29		
17.604.30	100	29	14	29	29		
17.604.35	100	45	14	14	14		
17.604.45	100	45	14	29	29		
17.605.20	100	29	14	29	29	29	
17.605.25	100	29	14	14	29	29	
17.605.30	100	45	14	29	29	29	
17.605.35	100	45	14	14	29	29	

Blank Lid

1-Neck

2-Neck

3-Neck


4-Neck

5-Neck

Special designs on request

Flat Flange Lid DN 120

Flat Flange SCHMIZO without a groove

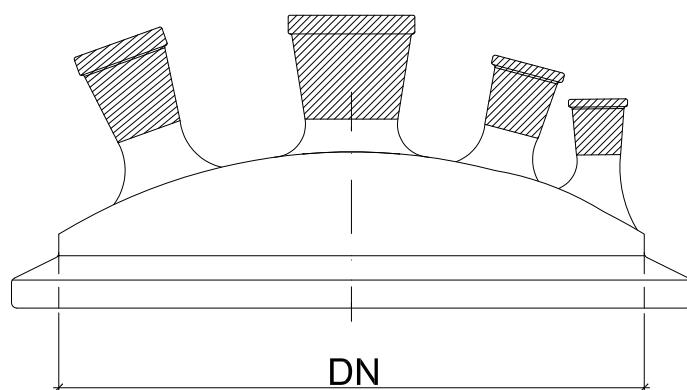

SCHMIZO Flat Flange

DN = Diameter nominal

All lids are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number.

Article No.	DN	Central NS	Lateral Vertical NS	Lateral Angled NS	Lateral Angled NS	Lateral Angled NS	
17.700.00	120						
17.701.10	120	14					
17.701.25	120	29					
17.701.35	120	45					
17.701.45	120	71					
17.702.10	120	14	14				
17.702.25	120	29	14				
17.702.30	120	29	29				
17.702.35	120	45	14				
17.702.40	120	45	29				
17.703.20	120	29	14	14			
17.703.25	120	29	14	29			
17.703.30	120	45	14	14			
17.703.40	120	45	14	29			
17.704.20	120	29	14	14	14		
17.704.25	120	29	14	14	29		
17.704.30	120	29	14	29	29		
17.704.35	120	45	14	14	14		
17.704.45	120	45	14	29	29		
17.705.20	120	29	14	29	29	29	
17.705.25	120	29	14	14	29	29	
17.705.30	120	45	14	29	29	29	
17.705.35	120	45	14	14	29	29	

Blank Lid

1-Neck

2-Neck

3-Neck


4-Neck

5-Neck

Special designs on request

Flat Flange Lid DN 150

Flat Flange SCHMIZO without a groove

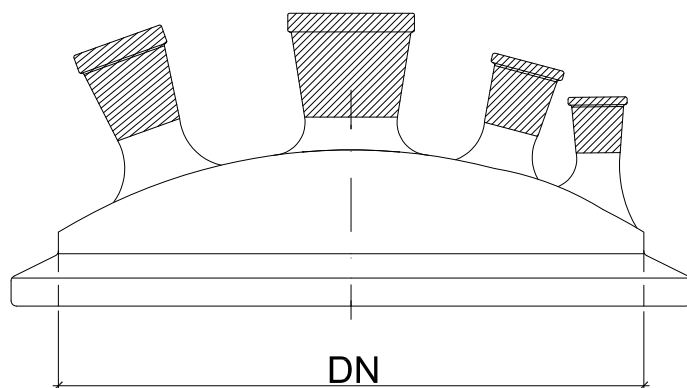

SCHMIZO Flat Flange

DN = Diameter nominal

All lids are also available with „SCHOTT“ flange (DURAN laboratory flange): Please add „S“ at the end of the article number.

Article No.	DN	Central NS	Lateral Vertical NS	Lateral Angled NS	Lateral Angled NS	Lateral Angled NS	
17.800.00	150						
17.801.10	150	14					
17.801.25	150	29					
17.801.35	150	45					
17.801.45	150	71					
17.802.10	150	14	14				
17.802.25	150	29	14				
17.802.30	150	29	29				
17.802.35	150	45	14				
17.802.40	150	45	29				
17.803.20	150	29	14	14			
17.803.25	150	29	14	29			
17.803.30	150	45	14	14			
17.803.40	150	45	14	29			
17.804.20	150	29	14	14	14		
17.804.25	150	29	14	14	29		
17.804.30	150	29	14	29	29		
17.804.35	150	45	14	14	14		
17.804.45	150	45	14	29	29		
17.805.20	150	29	14	29	29	29	
17.805.25	150	29	14	14	29	29	
17.805.30	150	45	14	29	29	29	
17.805.35	150	45	14	14	29	29	

Blank Lid

1-Neck

2-Neck

3-Neck


4-Neck

5-Neck

Special designs on request

Flat Flange Lid DN 200

Flat Flange SCHOTT without a groove

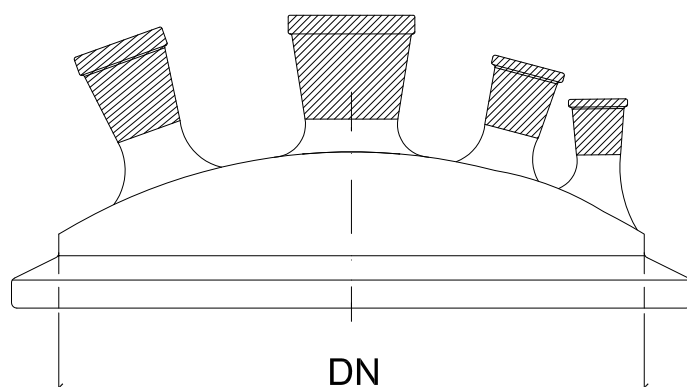

SCHMIZO Flat Flange

DN = Diameter nominal

All lids DN 200 are only available with „SCHOTT“ flange (DURAN laboratory flange).

Article No.	DN	Central NS	Lateral Vertical NS	Lateral Angled NS	Lateral Angled NS	Lateral Angled NS	
17.900.00S	200						
17.901.25S	200	29					
17.901.35S	200	45					
17.901.45S	200	71					
17.902.25S	200	29	14				
17.902.30S	200	29	29				
17.902.35S	200	45	14				
17.902.40S	200	45	29				
17.903.20S	200	29	14	14			
17.903.25S	200	29	14	29			
17.903.30S	200	45	14	14			
17.903.40S	200	45	14	29			
17.904.20S	200	29	14	14	14		
17.904.25S	200	29	14	14	29		
17.904.30S	200	29	14	29	29		
17.904.35S	200	45	14	14	14		
17.904.45S	200	45	14	29	29		
17.905.20S	200	29	14	29	29	29	
17.905.25S	200	29	14	14	29	29	
17.905.30S	200	45	14	29	29	29	
17.905.35S	200	45	14	14	29	29	

Blank Lid

1-Neck

2-Neck

3-Neck

4-Neck

5-Neck

Special designs on request

Filter Element & Filter Support

For filter funnel reactors, a filter element is clamped between a column extender and a glass part with outlet and a bottom valve. This setup allows to easily replace the filter whenever required.




Bottom Part (6), fits Filter Part and Filter Support with flat flange SCHMIZO without groove, glass part for bottom valve Schmizo *

Article No.	DN	Outlet Ø mm	
23.016.60	60	10	
23.016.61	100	10	
23.016.63	150	10	
23.016.64	200	10	
23.016.93	150	20	
23.016.94	200	20	

* also available with SCHOTT-flange: Please add „S“ at the end of the article number.

other valves on request

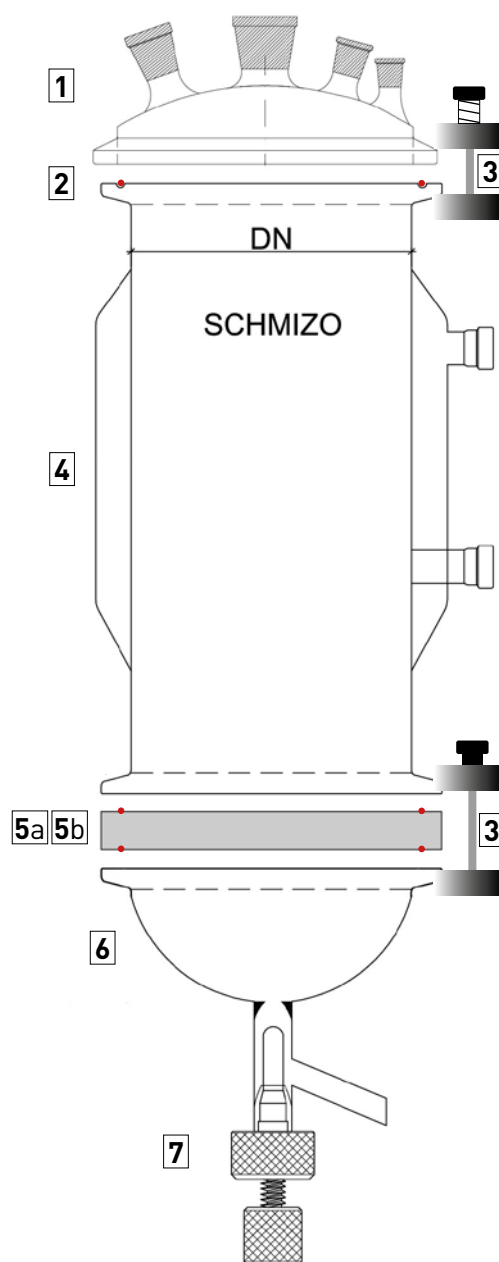

Design of a Filter Funnel Reactor:

- 1 Cover (from page 29)
- 2 O-Ring (s. 337)
- 3 2x reactor holder (s. 38)
- 4 Column extender (p. 26/27)
- 5 Filter part / Filter support
O-Rings included
- 6 Bottom part with outlet
- 7 Outlet valve (s. 40)

Filter Element PTFE (5a) with Sealed Glass Filter Plate

O-Rings are included

$x = \text{Porosity } [0 - 5]$

Article No. Filter element	DN	
23.016.06x	60	
23.016.10x	100	
23.016.15x	150	
23.016.20x	200	



Filter element

Filter Support PTFE (5b) for Paper or Membrane Filters

O-Rings are included

Article No. Filter support	DN	
23.016.50	60	
23.016.51	100	
23.016.53	150	
23.016.54	200	



Filter support

Filter Elements and Filter Support are not available for the size DN 120

Press Ring PTFE for Paper and Membran Filters

O-Rings are included

When using thicker filter papers or membranes, we recommend using a pressure ring to ensure tightness.

The pressure ring must be matched to the thickness of the filter in order to function optimally.

Please specify the thickness of the filter material used in your inquiry.



Filter support with filter and press ring

Reactor with Glass Filter Plate

All reaction vessels up to a volume of 10 liters max. - if single, jacketed or vacuum jacketed reactors - can also be manufactured with fixed glass filters of all porosities.

We would be pleased to provide an offer according to your requirements.



Jacketed reactor with fused-in glass filter plate P3


Volumes of Filter Reactors

The volume of filter reactors can be understood in two ways. Please specify the desired option you are interested in.


— Glass filter plate

Volume of the Reactor

Option 1

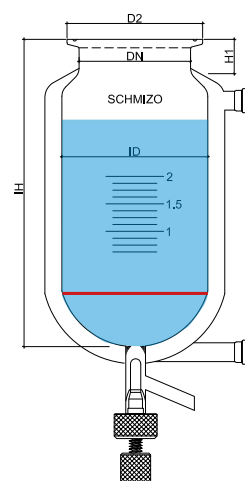

Option 2

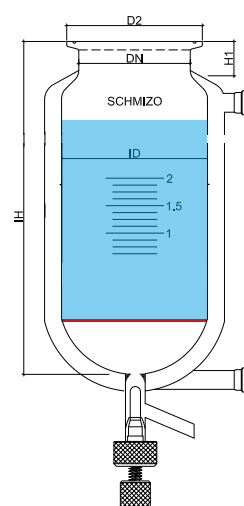

International Standards for Porosities in Glass Sintered Filters

ISO International Standard

ASTM American Standard

BS British Standard

We usually offer filter plates according to International Standard. Other porosities on request.

Porosities according to ISO 4793-80

Porosity	Pore size μm
Por. 00	250 - 500
Por. 0	160 - 250
Por. 1	100 - 160
Por. 2	40 - 100
Por. 3	16 - 40
Por. 4	10 - 16
Por. 5	1 - 1,6

Seals & Holders for Flat Flanges

O-Ring FEP covered with Silicone Core (-80°C up to 200°C)

Flat Flange **SCHMIZO**

Article No.	DN	
17.930.10	60	
17.930.20	100	
17.930.30	120	
17.930.40	150	

Flat Flange **SCHOTT**


Article No.	DN	
17.935.10	60	
17.935.20	100	
17.935.30	120	
17.935.40	150	
17.935.50	200	



O-Ring Viton/FPM 75 Shore (-20°C up to 200°C)

Flat Flange **SCHMIZO**

Article No.	DN	
17.920.10	60	
17.920.20	100	
17.920.30	120	
17.920.40	150	

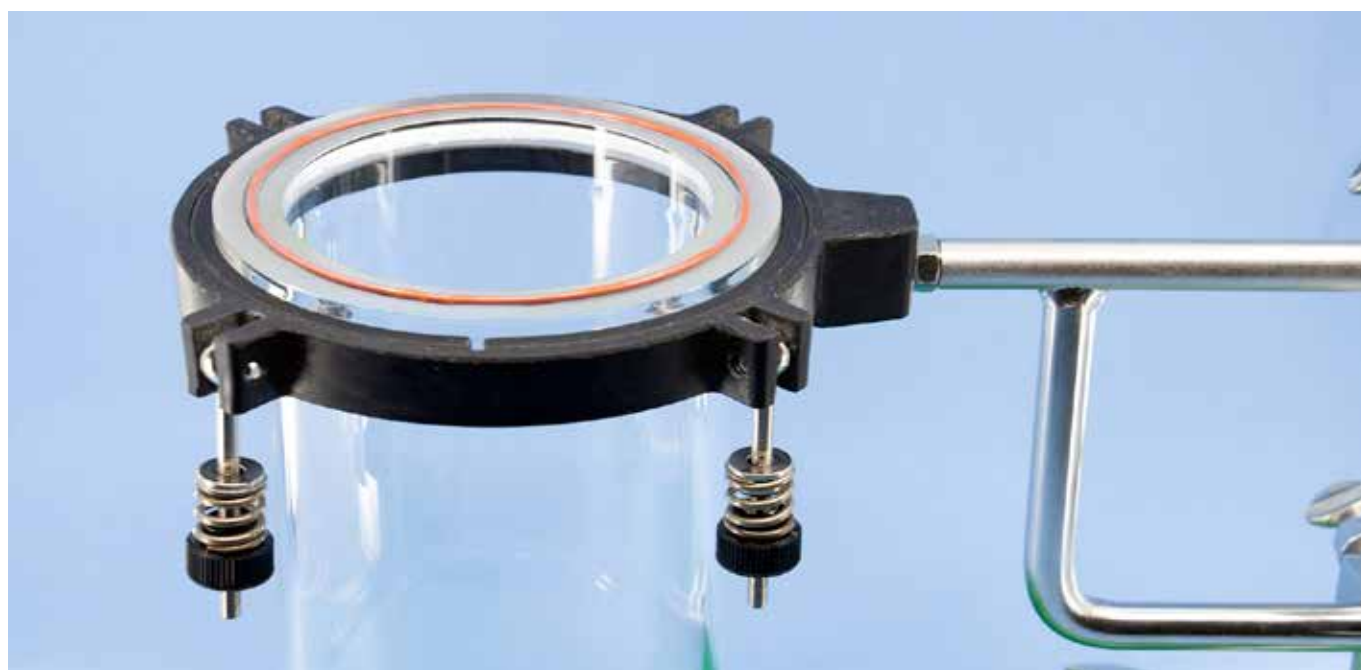

Quick Release Clamb, Stainless Steel

Article No.	DN	Fits Flange Type	
12.343.60	60	SCHMIZO / SCHOTT	
12.343.10	100	SCHMIZO / SCHOTT	
12.343.12	120	SCHOTT	
12.343.15	150	SCHMIZO / SCHOTT	
12.343.20	200	SCHOTT	



Reactor Holder Type SCHMIZO, Polyester (glass-fiber reinforced) , Max. Temp. 160° C (temporary up to 200° C)

Article No.	DN	Fits Flange Type	
12.344.98	60	SCHMIZO	
12.346.98	60	SCHOTT	
12.344.37	100	SCHMIZO	
12.346.37	100	SCHOTT	
12.344.50	120	SCHMIZO	
12.346.50	120	SCHOTT	
12.344.80	150	SCHMIZO	
12.346.80	150	SCHOTT	



Picture1: The reactor holders DN 60 / DN 100 / DN 120 are delivered with a thread rod M8 13 x 200 mm. Optionally you can modify them with a twin rod (pict. 3).

Picture1: The reactor holder DN 150 has two lateral support segments for rods up to 13 mm diameter. They are not included in delivery. See also picture following page.

Twin Rod for Reactor Holder Type SCHMIZO

Article No.	DN	Fits Flange Type	
12.348.40	60 - 120	SCHMIZO / SCHOTT	

Picture 3: The twin rod can be ordered optionally. It increases stability when carrying heavy loads. The horizontal alignment of the reactor holder is secured with a lock nut.

Reactor Holder, Stainless Steel

for lateral fixation with two rods $\varnothing$ 16 mm. The tension is transferred via the two spiral spring rings when the screws are tightened.

The support rods (see also page 68) are not included in the scope of delivery.

Article No.	DN	Fits Flange Type	
12.349.10	100	SCHOTT	
12.349.12	120	SCHOTT	
12.349.15	150	SCHOTT	
12.349.20	200	SCHOTT	

Picture 1:
Reactor holder stainless steel DN 100 with 4 support segments.

Picture 2:
Reactor holder stainless steel DN 200 with 6 support segments.



SCHMIZO reactor holder DN 150

Bottom Discharge Valve SCHMIZO

Bottom Discharge Valve SCHMIZO [-30°C up to +200°C]

Article No.	For Aperture mm	Construction	
17.940.10	10	standard	
17.940.12	10	O-Ring KALREZ	
17.940.14	10	O-Ring Viton	
17.940.16	10	with Pt-100 probe	
17.942.10	16	standard	
17.942.12	16	O-Ring KALREZ	
17.942.14	16	O-Ring Viton	
17.942.16	16	with Pt-100 probe	
17.944.10	20	standard	
17.944.12	20	O-Ring KALREZ	
17.944.14	20	O-Ring Viton	
17.944.16	20	with Pt-100 probe	

Also available with pneumatic cylinder.

The discharge valve can be used with limited leak tightness up to -100°C.



Construction of valve plug →



Standard

O-Ring

Pt-100

Components / Spare Parts

No.	Name
1	Valve plug PTFE
2	Fixing screw
3	Valve knob with thread and integrated pressure spring
4	Conical seal seat with O-ring joint
5	Spacer
6	Coupling nut
7	Spring

We can check if it is worth to revise a damaged discharge valve. The repair can be carried out in our workshop. The valve plug can only be replaced in our workshop.

The spare parts no. 2 - 7 can also be ordered separately. Information about the size of the plug will be necessary.

A manual with detailed instructions for assembly and disassembly of the bottom valve is available.



Hose Adapters for Jacket Inlet/Outlet

Stainless Steel Adapters for Flat Flanges DN 15 and DN 25
with PEEK connecting screw and O-Ring FEP/silicone

Article No.	Type	Connection	Form	Seal	
17.981.96P	DN 15	Olive 10mm	straight	O-Ring	
17.980.96P	DN 15	Olive 12mm	straight	O-Ring	
17.980.95P	DN 15	M 16x1	straight	O-Ring	
17.980.97P	DN 15	M 16x1	90°	O-Ring	
17.980.94P	DN 15	M 24x1.5	straight	O-Ring	
17.980.98P	DN 15	M 24x1.5	90°	O-Ring	
17.981.97P	DN 25	M 16x1	90°	O-Ring	
17.981.94P	DN 25	M 24x1.5	straight	O-Ring	
17.981.98P	DN 25	M 24x1.5	90°	O-Ring	

Type DN 25 with connection M30 and spare connection screws on request



Polyetheretherketon (abbreviated PEEK) is a high-temperature-resistant thermoplastic belonging to the polyaryletherketone group of substances. The melting temperature is 335 °C [Wikipedia].

Replacement Seals for Stainless Steel Adapter DN 15 / DN 25

Art.Nr.	Typ	Grösse	Material	max. Temp.	
17.982.15	O-Ring	DN 15	FEP/Silicone	205°C *	
17.982.25	O-Ring	DN 25	FEP/Silicone	205°C *	

* Kalrez O-ring for higher temperatures on request



The new version with PEEK screw no longer needs a safety spring. The sealing is still done by an O-ring. The metal part is identical to the version with the stainless steel screw.

Another advantage is that the PEEK screw can remain on the reactor when the metal adapter is separated from the reactor.




Stainless Steel Adapter for Spherical Joints KS 19 and KS 29
 with metal connecting screw and PTFE seal

Article No.	Type	Thread	Form	Seal	
17.980.10	KS 19	M 16x1	straight	KS-D	
17.980.11	KS 19	M 16x1	90°	KS-D	
17.980.20	KS 29	M 16x1	straight	KS-D	
17.980.21	KS 29	M 16x1	90°	KS-D	
17.980.22	KS 29	M 24x1.5	straight	KS-D	
17.980.23	KS 29	M 24x1.5	90°	KS-D	


Seals and Springs for Spherical Joints Stainless Steel Adapters

Article No.	Type	Size	Material	max. Temp.	
12.355.22	Seal	KS 19	PTFE	260°C	
12.355.40	Seal	KS 29	PTFE	260°C	
17.982.55	Spring	KS 19	Edelstahl	260°C	
17.982.56	Spring	KS 29	Edelstahl	260°C	

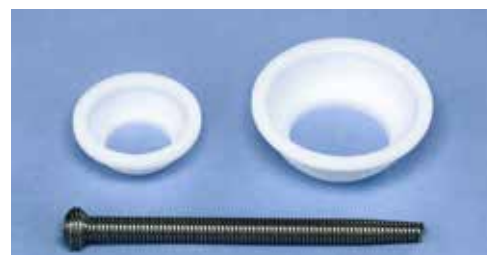
Spherical Joint with Hose Connection (borosilicate glass)

Article No.	SJ	
17.950.19	19/9	
17.950.29	29/15	

Other sizes are available on request

Clamps for Spherical Joints
 stainless steel, with tightening screw

Article No.	Size SJ	
12.341.14	19	
12.341.18	29	


PTFE Hose Connection with PBT Cap for GL-Threads

Article No.	Form	for GL	ID mm	OD mm	
12.259.40	straight	14	5.0	9.0	
12.259.42	straight	18	7.5	10.5	
12.259.44	straight	25	10	13	
12.259.46	arced	14	5.0	9.0	
12.259.48	arced	18	7.5	10.5	



PP Hose Connection with PBT Cap for GL-Threads

Article No.	Form	for GL	ID mm	OD mm	
12.259.03	straight	14	5.0	8.6	
12.259.06	straight	18	6.8	10.8	
12.259.05	arced	14	4.0	8.6	
12.259.08	arced	18	6.8	10.8	



Insulated Metal Hoses

Metal Hose STANDARD

Wave hose made of stainless steel with special silicone foam insulation and end cap, fitting SCHMIZO reactors and hose adapters.

Article No.	Thread	Length m	DN mm	Temperature range °C	
17.985.10	M 16x1	1.0	10	-100 ... +260	
17.985.15	M 16x1	1.5	10	-100 ... +260	
17.985.20	M 16x1	2.0	10	-100 ... +260	
17.985.30	M 16x1	3.0	10	-100 ... +260	
17.985.55	M 24x1.5	1.0	12	-100 ... +260	
17.985.60	M 24x1.5	1.5	12	-100 ... +260	
17.985.65	M 24x1.5	2.0	12	-100 ... +260	
17.985.70	M 24x1.5	3.0	12	-100 ... +260	



Metal Hose ECO

Wave hose made of stainless steel with silicone insulation, fitting SCHMIZO reactors and hose adapters.

Article No.	Thread	Length m	DN mm	Temperature Range °C	
17.986.10	M 16x1	1.0	10	-30 ... +200	
17.986.15	M 16x1	1.5	10	-30 ... +200	
17.986.20	M 16x1	2.0	10	-30 ... +200	
17.986.30	M 16x1	3.0	10	-30 ... +200	




Anchor Stirrers

KPG Anchor-Stirrer (KPG-length 160 mm)

Please indicate desired length **L** when placing your order!

Article No.	Ø KPG mm	Chuckling Ø mm	Width B mm	Height H mm	
24.603.06	10	8	60	60	
24.603.09	10	8	90	90	
24.603.11	10	8	110	110	
24.603.14	10	8	140	140	
24.604.06	16	10	60	60	
24.604.09	16	10	90	90	
24.604.11	16	10	110	110	
24.604.14	16	10	140	140	


*All stirrers are also available with ZRS instead of KPG
(ZRS-length 300mm).*

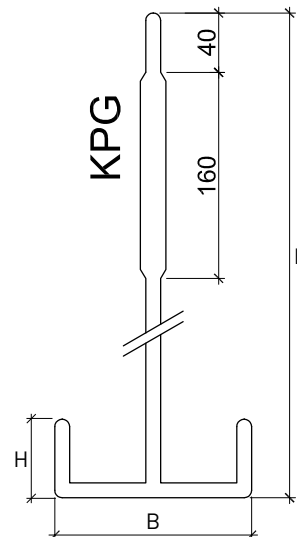

Anchor Stirrer with Simple Glass Shaft

Please indicate desired length **L** when placing your order!

Article No.	Shaft Ø mm	Chuckling Ø mm	Width B mm	Height H mm	
24.605.06	8	8	60	60	
24.605.09	8	8	90	90	
24.605.11	8	8	110	110	
24.605.14	8	8	140	140	
24.606.06	10	10	60	60	
24.606.09	10	10	90	90	
24.606.11	10	10	110	110	
24.606.14	10	10	140	140	

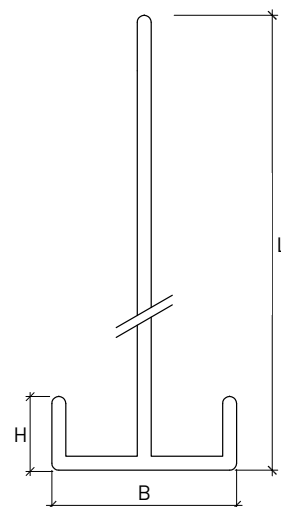

auch als Doppel-Anker erhältlich


4 Blade Stirrer

KPG 4 Blade Stirrer (KPG Length 160 mm)

Please indicate desired length **L** when placing your order!

Article No.	Ø KPG mm	Chucking Ø mm	Width B mm	Height H mm	
24.634.24	10	8	24	25	
24.634.40	10	8	40	25	
24.634.60	10	8	60	25	
24.635.24	16	10	24	25	
24.635.40	16	10	40	25	
24.635.60	16	10	60	25	


All stirrers are also available with ZRS instead of KGP (ZRS-length 300mm).

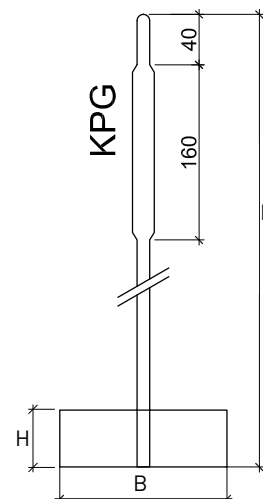

4 Blade Stirrer with Simple Glass Shaft

Please indicate desired length **L** when placing your order!

Article No.	Shaft Ø mm	Chucking Ø mm	Width B mm	Height H mm	
24.636.24	8	8	24	25	
24.636.40	8	8	40	25	
24.636.60	8	8	60	25	
24.637.24	10	10	24	25	
24.637.40	10	10	40	25	
24.637.60	10	10	60	25	

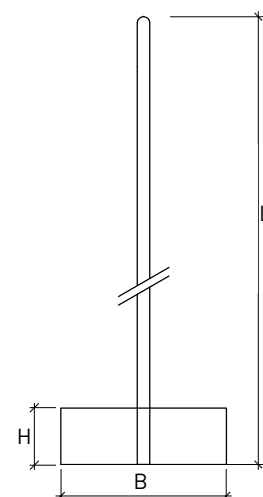

auch als Doppel-Version erhältlich

Stirrer Shafts with Eye


for PTFE blades with bolt and straddling stirrer

KPG Stirrer Shaft with Eye (KPG-Length 160 mm)

Please indicate desired length **L** when placing your order!

Article No.	Ø KPG mm	Chuckling Ø mm	
23.638.21	10	8	
23.638.24	16	10	


All stirrers are also available with ZRS instead of KGP (ZRS-length 300mm).

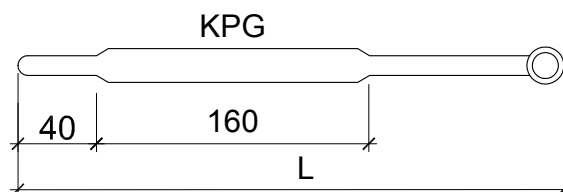

Stirrer shaft with Eye with Simple Glass Shaft

Please indicate desired length **L** when placing your order!

Article No.	Ø shaft mm	Chuckling Ø mm	
23.629.08	8	8	
23.629.10	10	10	
23.629.12	12	10	
23.629.14	14	10	
23.629.16	16	10	

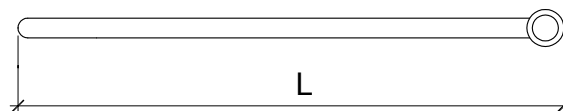


Straddling Stirrer Blade with Bolt

for stirrer shafts with eye, suitable for joints size NS 24 or bigger.

Article No.	Stirring Ø mm	Bolt Ø mm	
23.643.07	50	8	
23.643.09	70	8	
23.643.11	90	8	

*Replacement bolts and other sizes on request.
(Available with blade widths 30, 120 and 140 mm,
however with different bolt sizes!)*

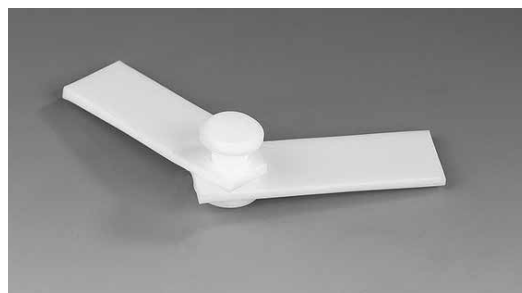

Stirrer Shafts with Double Pivot


for PTFE moon-shaped stirrer blades with double-sided eyelet

KPG Stirrer Blade with Double Pivot (KPG Length 160 mm)

Please indicate desired length **L** when placing your order!

Article No.	Ø KPG mm	Chuckling Ø mm	
23.638.01	10	8	
23.638.04	16	10	


All stirrers are also available with ZRS instead of KGP (ZRS-length 300mm).

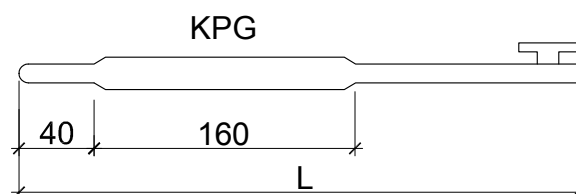

Stirrer Blade with Double Pivot with Simple Glass Shaft

Please indicate desired length **L** when placing your order!

Article No.	Ø Shaft mm	Chuckling Ø mm	
23.640.01	10	10	

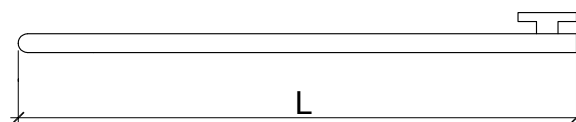

PTFE Moon-Shaped Stirrer Blades with Double-Sided Eyelet

for stirrer shafts with double pivot

Article No.	Stirrer blade mm	fits NS	
23.644.01	50 x 18 x 3	24	
23.644.02	50 x 24 x 3	29	
23.644.04	65 x 18 x 3	24	
23.644.05	68 x 24 x 3	29	
23.644.06	75 x 24 x 3	29	
23.644.07	90 x 24 x 3	29	
23.644.08	110 x 24 x 3	29	
23.644.09	125 x 24 x 3	29	
23.644.10	145 x 35 x 4	45	



Stirrer Shaft with Thread Head


for PTFE stirrer heads with SVL TORION 5 thread

KPG Stirrer Shaft with Thread Head (KPG Length 160 mm)

Please indicate desired length **L** by placing your order!

Article No.	Ø KPG mm	Chuckling Ø mm	
23.647.01	10	8	
23.647.04	16	10	


All stirrers are also available with ZRS instead of KGP (ZRS-length 300mm).

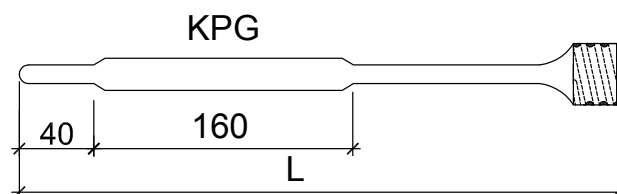

Stirrer Shaft with Thread Head with Simple Glass Shaft

Please indicate desired length **L** by placing your order!

Article No.	Ø Shaft mm	Chuckling Ø mm	
23.648.01	10	10	
23.648.02	12	10	
23.648.03	14	10	
23.648.04	16	10	

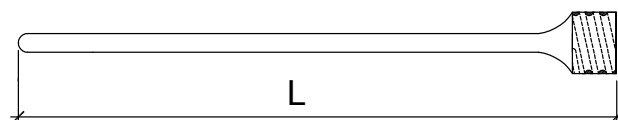

PTFE Stirrer Heads with Thread

fits stirrer shafts with SVL TORION 5 thread

Article No.	Stirrer head type	Wide mm	
23.650.10	Straddling blades (1)	70	
23.650.20	Leaf (2)	70	
23.650.30	Small anchor (3)	54	
23.650.35	Big anchor (4)	90	



1



2



3




4

Stainless Steel Stirrer


Anchor Stirrer

Article No.	Width B mm	Height H mm	Length L mm	Ø Shaft mm	
23.651.01	60	60	800	10	
23.651.02	60	60	1000	10	
23.651.41	90	90	800	10	
23.651.42	90	90	1000	10	
23.651.03	110	110	800	10	
23.651.04	110	110	1000	10	
23.651.05	140	140	800	10	
23.651.06	140	140	1000	10	
23.651.07	140	140	1200	10	
23.651.08	190	190	800	10	
23.651.09	190	190	1000	10	
23.651.10	190	190	1200	10	

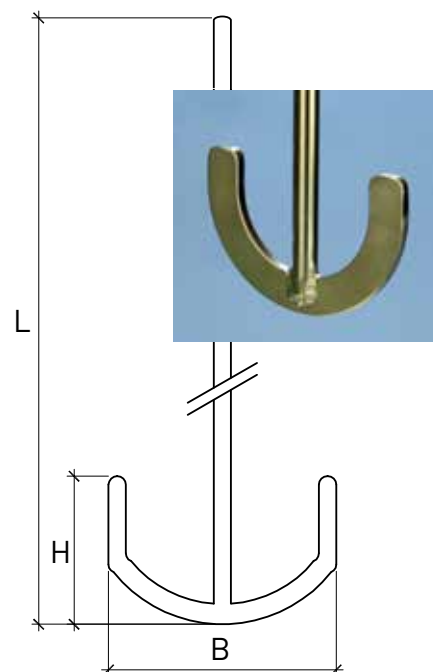

Pretzel-Stirrer

Article No.	Width B mm	Height H mm	Length L mm	Ø Shaft mm	
23.651.21	60	60	800	10	
23.651.22	60	60	1000	10	
23.651.45	90	90	800	10	
23.651.46	90	90	1000	10	
23.651.23	110	110	800	10	
23.651.24	110	110	1000	10	
23.651.25	140	140	800	10	
23.651.26	140	140	1000	10	
23.651.27	140	140	1200	10	
23.651.28	190	190	800	10	
23.651.29	190	190	1000	10	
23.651.30	190	190	1200	10	

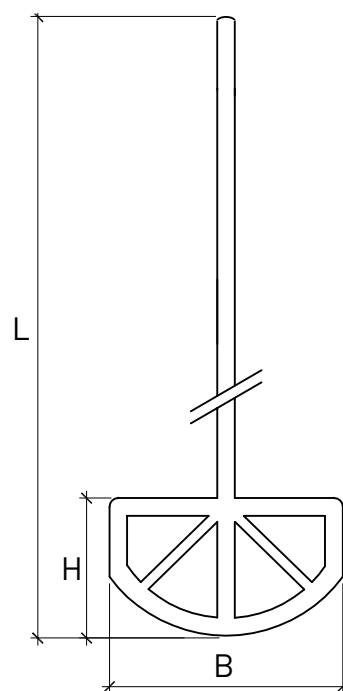

All other forms and sizes by request. Stainless steel stirrer shafts can be coated with HALAR (max. working temperature- 150°C).

BOLA Stirrers Stainless Steel/PTFE

As alternative to glass or stainless steel stirrers, we offer BOLA stirrers. Please find the whole BOLA range on:


<http://www.bola.de/en/products/stirring-and-mixing/stirrer-shafts.html>

BOLA stirrers are made with a PTFE covered stainless steel shaft and blades of PTFE. The steel core provides high stability and allows safe fixing in the stirrer coupling. The product which is stirred is only exposed to PTFE.

Temperature range: -200°C up to +250°C

Chemical resistance: +++ universal

Example: **Anchor Stirrer**

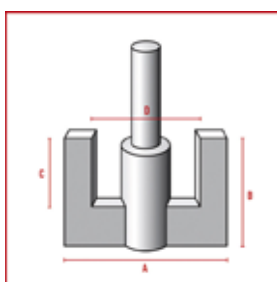

BOLA Stirrer with one Paddle



BOLA Stirrer with two Paddles



BOLA Propeller Stirrer



BOLA Centrifugal Stirrer



BOLA Disc Stirrer

Article No.	Length mm	Shaft Ø mm	Chuckling Ø mm	A mm	B mm	C mm	D mm
C 384-01	350	8	6.5	40	35	20	26
C 384-02	350	8	6.5	60	40	25	30
C 384-04	450	8	6.5	60	40	25	30
C 384-06	450	8	6.5	80	50	30	44
C 384-08	600	8	6.5	80	50	30	44
C 384-10	600	8	6.5	100	60	35	56
C 384-16	350	10	8	80	50	30	44
C 384-17	450	10	8	80	50	30	44
C 384-24	600	10	8	100	60	35	56
C 384-28	800	10	8	100	60	35	56
C 384-32	1000	10	8	100	60	35	56
C 384-40	1200	10	8	100	60	35	56
C 384-44	600	10	8	130	80	55	80
C 384-48	800	16	8	130	80	55	80
C 584-52	800	16	10	150	120	90	90
C 584-58	1000	16	10	150	120	90	90
C 584-64	1200	16	10	150	120	90	90

Stirrer Bearings with Glass Shaft Guide

KPG Stirrer Bearing with NS Joint Cone, optionally with Grease Cup
suitable for KPG stirrer bearings

Article No.	Ø KPG mm	NS Joint Cone	grease cup	
23.601.32	10	24	-	
23.602.32	10	29	-	
23.603.32	10	45	-	
23.604.32	10	24	✓	
23.605.32	10	29	✓	
23.606.32	10	45	✓	
23.606.48	16	45/40	✓	



Jacketed Stirrer Bearing

suitable for simple glass or KPG stirrer shafts

The screw cap with seal for the upper thread and the screw caps with hoses for the Inlet/Outlet are included.

Article No.	Inner Ø	NS Joint Cone	Upper thread	Inlet/ Outlet	
23.610.32	8 mm	24	GL 18	GL 14	
23.610.40	8 mm	29	GL 18	GL 14	
23.610.42	8 mm	45	GL 18	GL 14	
23.610.46	8 mm	71	GL 18	GL 14	
23.612.32	KPG 10	24	GL 18	GL 14	
23.612.40	KPG 10	29	GL 18	GL 14	
23.612.48	KPG 10	45	GL 18	GL 14	
23.615.45	KPG 16	45	GL 25	GL 14	



Jacketed KPG Stirrer Bearing, KPG Inside the Jacket

suitable for KPG stirrer shafts.

Screw caps with hoses are included.

Article No.	KPG Ø mm	NS Joint Cone	Inlet/ Outlet	
23.613.01	10	29	GL 14	
23.613.02	10	45	GL 14	
23.613.23	16	45	GL 14	
23.613.24	16	60	GL 14	
23.613.25	16	71	GL 14	



Jacketed KPG Stirrer Bearing, KPG above the Jacket

suitable for KPG stirrer shafts.

Srew caps with hoses are included.

Article No.	Ø KPG mm	NS Joint Cone	Inlet/ Outlet	
23.613.11	10	29	GL 14	
23.613.12	10	45	GL 14	
23.613.33	16	45	GL 14	
23.613.44	16	60	GL 14	
23.613.45	16	71	GL 14	



Stirrer Bearing Eco

Stirrer Bearing with PTFE Shaft Guide

suitable for all kind of stirrer shafts. Gasket and guide shaft are exchangeable.

Article No.	Ø mm	NS Joint Cone	Upper Thread	
23.616.10	8	29	GL 25	
23.616.15	10	29	GL 25	
23.616.20	10	45	GL 25	
23.616.25	16	45	GL 32	



Replacement PTFE Shaft Guide, Gasket included

Article No.	Ø mm	
23.616.01	8	
23.616.02	10	
23.616.03	16	



Replacement Special Gasket PTFE/FPM

Article No.	Ø mm	
23.616.11	8	
23.616.12	10	
23.616.13	16	



Stirrer Bearing DURA

Vacuum* and pressure tight due to variable, easily replaceable sealing elements made of PTFE (brown parts with glass fiber additive).


* Leaking rate at 400 rpm: $1,2 \times 10^{-3}$ mbar/s
(with new seal kits)

Stirrer Bearing DURA


with replaceable double lip seal

Article No.	Ø mm	NS Cone	
23.622.03	10	29	
23.622.04	10	45	
23.622.07	16	29	
23.622.08	16	45	

Spare Parts for Stirrer Bearing DURA

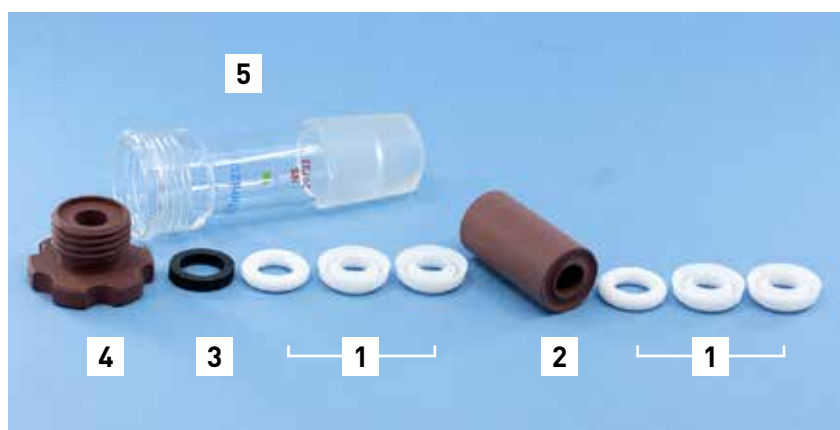

DURA is suitable for vacuum, pressure or universal use, depending on the installation of the sealing elements:

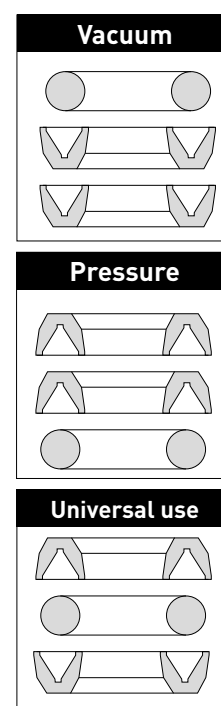

Article No.	Ø mm	Spare part	
23.622.10	10	Seal kit (1) *	
23.622.11	10	Guide sleeve (2)	
23.622.12	10	Square ring Viton (3)	
23.622.13	10	Pressure screw (4) with square ring	
23.622.17	10	Glass body NS 29/32 (5)	
23.622.18	10	Glass body NS 45/40 (5)	
23.622.20	16	Seal kit (1) *	
23.622.21	16	Guide sleeve (2)	
23.622.22	16	Square ring Viton (3)	
23.622.23	16	Pressure screw (4) with square ring	
23.622.27	16	Glass body NS 29/32 (5)	
23.622.28	16	Glass body NS 45/40 (5)	


* The complete variant contains 2 seal kits (1), consisting of 2 lip seals and 1 PTFE O-ring.

Stirrer Bearing FLUKA

Suitable for vacuum and lightly over-pressure

The core part of the FLUKA stirrer bearing is the shaft guide made of PTFE, which is equipped with a shaft seal and two inside FPM O-rings for guiding the stirrer rod. The embedded O-ring seals the guide sleeve against the glass part. In addition to the optimized guidance of the stirring rod, this guiding system enables to work under vacuum and slight overpressure.

The stirrer bearing is available either with a screw cap or a suction head.

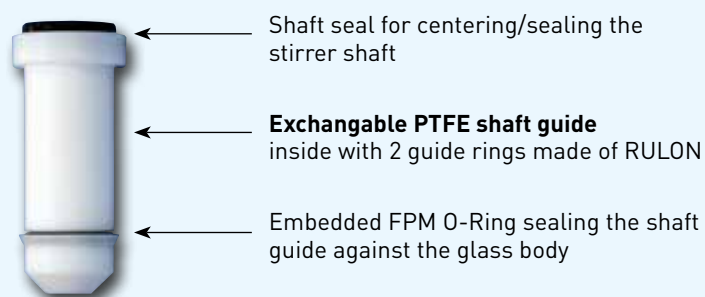

Stirrer Bearing FLUKA with GL-Screw Cap

with exchangeable PTFE shaft guide and shaft seal for centering. Suitable for all kind of stirrer shafts.

Article No.	Ø mm	NS Joint Cone	GL	
23.619.20	10	29	32	
23.619.22	10	45	32	
23.619.40	16	29	45	
23.619.42	16	45	45	

Due to lack of demand, we no longer offer the Ø 8mm and Ø 14mm variants.



Replacement Shaft Guide FLUKA

shaft seal included; suitable for all kind of stirrer shafts.

Article No.	Ø mm	
23.619.25	10	
23.619.45	16	



Replacement Shaft Seal Typ A, NBR

Quality according to DIN 3760 / ISO 6194

Article No.	Ø mm	
23.689.15	10	
23.689.35	16	



Stirrer Bearing FLUKA with Suction Head

with exchangeable PTFE shaft guide, centering with 2 shaft seals.
Fitting all kind of stirrer shafts.

Article No.	Ø mm	NS Joint Cone	GL	
23.618.20	10	29	32	
23.618.22	10	45	32	
23.618.40	16	29	45	
23.618.42	16	45	45	

Suction Head for Stirrer Bearing FLUKA

Article No.	GL	
23.618.00	32	
23.618.02	45	



Büchi Magnetic Couplings

For the highest demands, we offer two versions of Büchi magnetic couplings. The drive is contactless via a magnetic field, this is guaranteeing no abrasion at all. The seal to the glass cone is provided by two O-rings (FKM or FFKM).

The "bmd 075" models are driven by any external stirrer motor, while the "Cyclone 075 dc" models have an integrated drive.

With a maximum torque of 75 Ncm, these couplings are designed for operation under extreme conditions. The temperature range is from -50 °C to +200 °C. The magnetic couplings can withstand a full vacuum of -1 bar.

Perfect for vessels up to 5 liters.

Art.no.	Articel	
26.653.00	Magnetic Coupling bmd 075 stainless steel	
26.653.01	Magnetic Coupling bmd 075 Hastelloy®	
26.653.10	Magnetic Coupling Cyclone 075 dc stainless steel	
26.653.11	Magnetic Coupling Cyclone 075 dc Hastelloy®	

All models with NSK 29/32, stirrer holder Ø 10 mm



Mechanical Seals


Mechanical seals ensure a maximum of tightness with a minimum of abrasion.

Characteristics / Benefits:

- Completely oil and grease free
- Vacuum up to 1 mbar
- Excellent chemical resistance
- FDA compatible
- Meets GMP-guidelines
- Joint with tempering optionally
- Stainless steel ball bearing
- Provides non shedding of material into the reaction

Materials in contact with the product:

Borosilicate glass, Perfluororelastomer-FFKM, PTFE, PFA

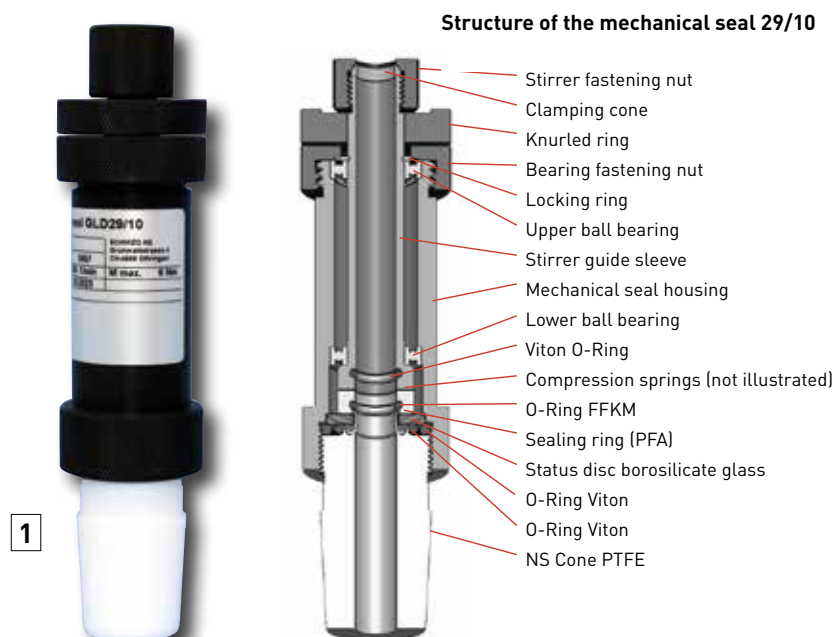

The models for stirrer shafts with $\varnothing 16$ mm, the NS glass cone is fixed with a spring and is equipped with a PTFE sleeve.

Spare parts, models for KF 25/3 and KF 40/3 as well as variants for the expulsion protected area on request.

Article No.	Ø Shaft mm	Cone	Installation Length mm	Material NS Cone	Pict no.	
23.623.16	8	NS 14	75	PTFE		
23.623.18	8	NS 19	75	PTFE		
23.623.20	10	NS 24	120	PTFE		
23.623.22	10	NS 29	120	PTFE	1	
23.623.24	10	NS 45	120	PTFE		
23.623.40	16	NS 45	160	Glass	2	
23.623.42 *	16	NS 45	160	Glass		

* Model with NS cone cooling

Stirrer Couplings

Stirrer Coupling with Double Cardan Joint

with compression fittings made of PVDF, handling without any tool. The PVDF head is exchangeable. Max. 1000 rpm.

Article No.	Opening Ø mm	Total length mm	
23.656.08	8	150	
23.656.10	10	150	



Replacement PVDF Head

Article No.	Opening Ø mm	
23.656.09	8	
23.656.11	10	



Clamp Coupling with Double Cardan Joint

especially for stainless steel stirrer shafts. Fixation of the stirrer shaft with an Allen screw. Max. 1000 rpm.

Article No.	Opening Ø mm	Total length mm	
23.656.18	8	130	
23.656.20	10	132	
23.656.21	14	134	
23.656.22	16	136	



BOLA Globe Stirrer Coupling made of POM plastics with compression fittings, handling without any tool. By turning the head the opening changes. Max. 1200 rpm.

Article No.	Opening Ø mm	Total length mm	
23.656.45	8 + 10	190	



CAT Stirrer Motors



CAT stirrer motors stand out with an excellent price-performance ratio.

Complete technical specifications you can find here: <http://www.cat-ing.de>



CAT R100 Series - high torque stirrer motors in four different executions for stirring volumes up to 200 liters of water.


CAT R20 - lightweight, maintenance-free stirrer motor for volumes up to 30 liters of water.

Article No.	Type	Engine output	Revolutions Rpm	Torque Ncm	Chuck mm	Weight kg	
23.658.10	R20	60 W	50 - 2000	20	0.5 - 8	1.3	
23.658.20	R100 SD	300 W	40 - 2000	100	1 - 10,5	5.0	
23.658.30	R100 C	300 W	40 - 2000	150	1 - 10,5	5.0	
23.658.40	R100 CL	300 W	30 - 1100	250	1 - 10,5	5.0	
23.658.50	R100 CT	300 W	10 - 500	600	1 - 10,5	5.5	

We also can offer stirrer motors from other manufacturers.

Liquid Dividers

The manually operated column heads are used to adjust the reflux ratio between reflux and drain in distillation columns. Depending on the version, they are designed for operation under vacuum and/or normal pressure.

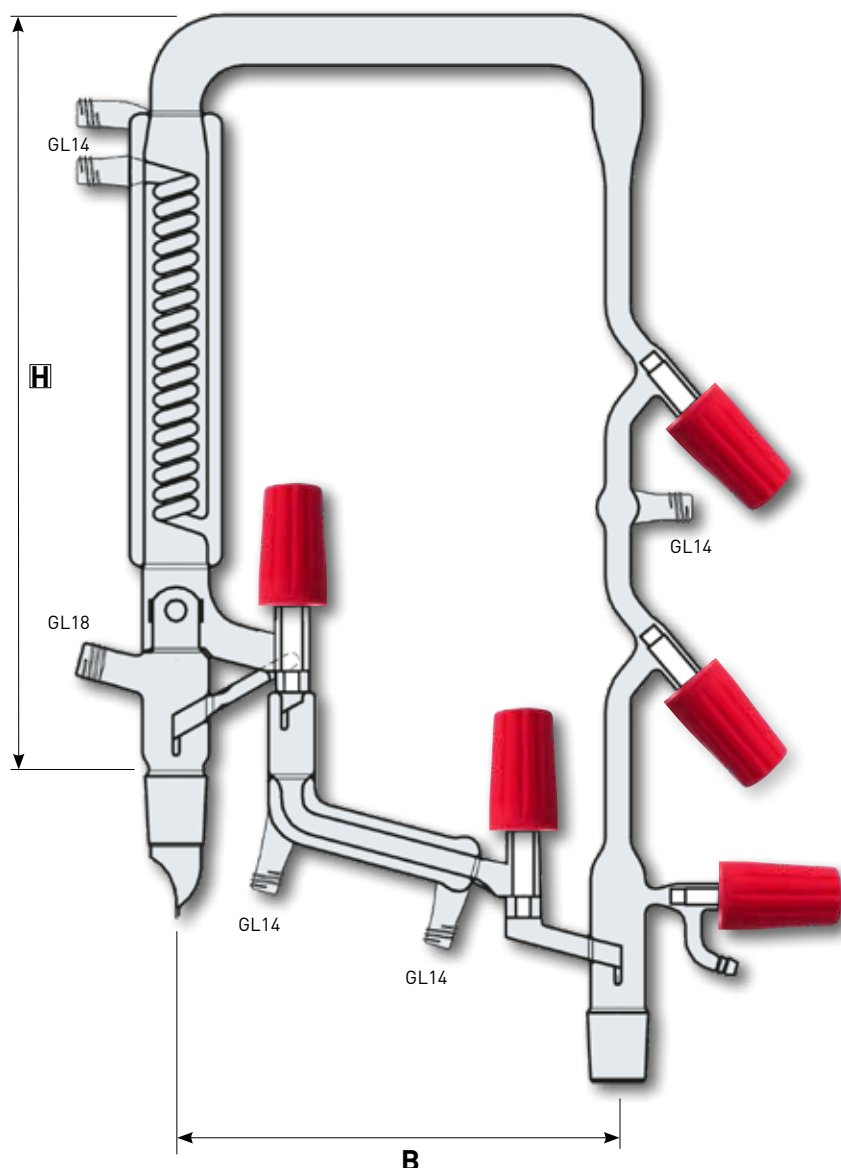

As there are countless variants and variations, we present here our best-selling liquid divider. Of course, all parameters of this model can also be configured according to your wishes. We will gladly offer you a liquid divider according to your needs.

PRODURAN® Valves



- Precise dosing
- Unscrew protection
- Vacuum-tight, leak rate $< 1 \times 10^{-6}$ mbar x l/s

Liquid Divider with Intensive and Intercooler, Dead Volume free

In- an outlet NS 29/32, intensive cooler 200 mm

Article No.	H mm	B mm	Valves	Vakuum operation	
20.738.30	370	170	Produran 6	yes	




All GL-14 connections with cap and PP olive (optionally PTFE)

Gas Inlet Tubes & Vacuum Pieces


Gas Inlet Tubes, right angle bend, with hose nozzle and NS joint.
Available with or without stopcock (glass/PTFE).
Please indicate desired tube length (as drawn).

Article No.	NS Joint Cone	Stopcock	
20.754.12	14/23	without	
20.754.40	29/32	without	
20.760.12	14/23	glass key, NS 14 / Bore 2.5	
20.760.40	29/32	glass key, NS 14 / Bore 2.5	
20.765.12	14/23	PTFE key, NS 14 / Bore 2.5	
20.765.40	29/32	PTFE key, NS 14 / Bore 2.5	

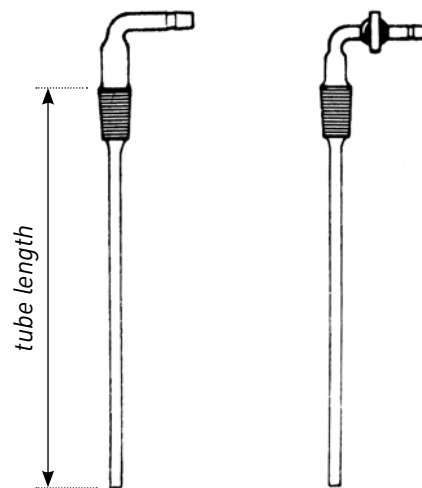


Vacuum Pieces Right Angle Bend, with hose nozzle and NS joint.
Available with or without stopcock (glass/PTFE).

Article No.	NS Joint Cone	Stopcock	
20.332.12	14/23	without	
20.332.40	29/32	without	
20.344.12	14/23	glass key, NS 14 / Bore 2.5	
20.344.40	29/32	glass key, NS 14 / Bore 2.5	
20.345.12	14/23	PTFE key, NS 14 / Bore 2.5	
20.345.40	29/32	PTFE key, NS 14 / Bore 2.5	
20.351.40	29/32	Produran valve 0 - 4 mm	

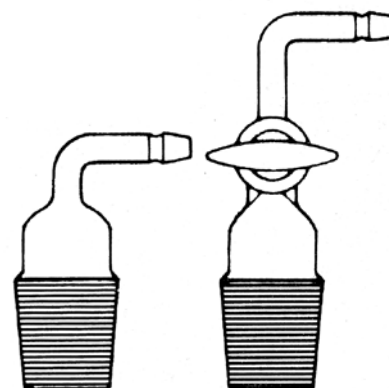

Vacuum Pieces Straight, with hose nozzle and NS joint.
Available with or without stopcock (glass/PTFE).


Article No.	NS Joint Cone	Stopcock	
20.330.12	14/23	without	
20.330.40	29/32	without	
20.342.12	14/23	glass key, NS 14 / Bore 2.5	
20.342.40	29/32	glass key, NS 14 / Bore 2.5	
20.343.12	14/23	PTFE key, NS 14 / Bore 2.5	
20.343.40	29/32	PTFE key, NS 14 / Bore 2.5	

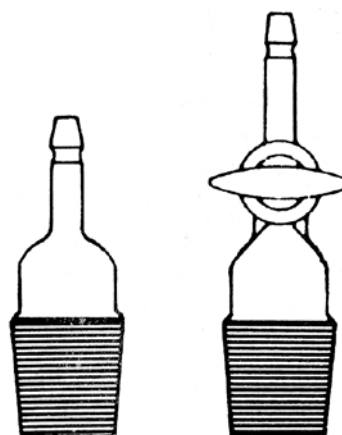

Vacuum Piece with High Vacuum Stopcock

laterally hose nozzle

Article No.	NS Joint Cone	Stopcock	
20.350.12	14/23	glass key, borehole 6 mm	
20.350.40	29/32	glass key, borehole 6 mm	

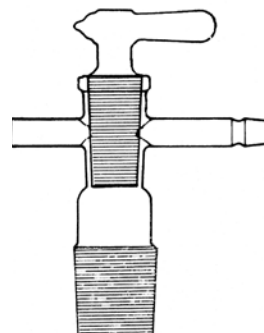


Alle articles are also available with NS-joint 19/26 or 24/29

Thermometer

Thermometer with Straight Stem, Stem Ø 8 mm, Subdivision 1°C

Article No.	Length mm	Measuring range	Filling	
40.100.15	150	-30 +150°C	Petroleum	
40.100.24	240	-30 +150°C	Petroleum	
40.100.33	330	-30 +150°C	Petroleum	
40.100.43	430	-30 +150°C	Petroleum	
40.100.60	600	-30 +150°C	Petroleum	
40.110.15	150	-30 +250°C	Isoamyl	
40.110.24	240	-30 +250°C	Isoamyl	
40.110.33	330	-30 +250°C	Isoamyl	
40.110.43	430	-30 +250°C	Isoamyl	
40.110.60	600	-30 +250°C	Isoamyl	
40.130.15	150	-100 +100°C	Toluol	
40.130.24	240	-100 +100°C	Toluol	
40.130.33	330	-100 +100°C	Toluol	
40.130.43	430	-100 +100°C	Toluol	
40.130.60	600	-100 +100°C	Toluol	

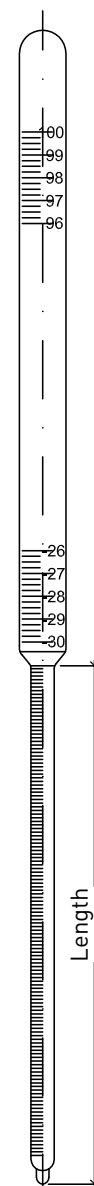

Temperature Probes

BOLA Temperature Probe - Lemo-Compact, PTFE

Thermocouple (PT 100) in a PTFE covered stainless steel tube.

Probe Ø 8 mm, tip Ø 6 mm, collar ring Ø 12 mm. Connection by a coupling (type Lemo®, socket size 1, 4-pole) fixed directly at the end of the probe.

Measuring range: -50°C to +250°C (Resistance -200°C to +250°C)



Article No.	Model	Total Length	Usable Length	
40.419.10	P 1730-10	150 mm	80 mm	
40.419.20	P 1730-20	370 mm	300 mm	
40.419.23	P 1730-23	470 mm	400 mm	
40.419.25	P 1730-25	570 mm	500 mm	

Probes made of stainless steel by request

Temperature Probe Pt-100 (Glass), acc. to DIN EN 60751, K1.A.

Tube Ø 8 mm made of borosilicate glass, head with screw cap with integrated Lemo® socket. Probe and tube are exchangeable.

Measuring range from -200°C to +400°C.



Article No.	Probe	Stock Length	Thread Head	
40.430.05	1 x Pt-100 2-Wire	150 mm	SVL 15	
40.430.10	1 x Pt-100 2-Wire	240 mm	SVL 15	
40.430.15	1 x Pt-100 2-Wire	330 mm	SVL 15	
40.430.20	1 x Pt-100 2-Wire	430 mm	SVL 15	
40.430.22	1 x Pt-100 2-Wire	600 mm	SVL 15	
40.430.45	1 x Pt-100 4-Wire	150 mm	GL 18	
40.430.50	1 x Pt-100 4-Wire	240 mm	GL 18	
40.430.55	1 x Pt-100 4-Wire	330 mm	GL 18	
40.430.60	1 x Pt-100 4-Wire	430 mm	GL 18	
40.430.62	1 x Pt-100 4-Wire	600 mm	GL 18	
40.430.65	2 x Pt-100 2-Wire	150 mm	GL 18	
40.430.70	2 x Pt-100 2-Wire	240 mm	GL 18	
40.430.75	2 x Pt-100 2-Wire	330 mm	GL 18	
40.430.80	2 x Pt-100 2-Wire	430 mm	GL 18	
40.430.82	2 x Pt-100 2-Wire	600 mm	GL 18	



SVL 15
Lemo 2-Pol

GL 18
Lemo 4-Pol

Replacement Tubes for Pt-100 Probes

made of borosilicate glass. Please indicate desired Probe Type

Article No.	Thread Head	Stock Length	
40.435.15	SVL 15	150 mm	
40.435.24	SVL 15	240 mm	
40.435.33	SVL 15	330 mm	
40.435.43	SVL 15	430 mm	
40.430.60	SVL 15	600 mm	
40.438.15	GL 18	150 mm	
40.438.24	GL 18	240 mm	
40.438.33	GL 18	330 mm	
40.438.43	GL 18	430 mm	
40.438.60	GL 18	600 mm	

black screw cap = SVL 15 (plug 2-pole)
red screw cap = GL 18 (plug 4-pole)



Guide tubes with other head threads by request

Replacement Probes Pt-100 According to DIN

incl. Lemo socket and screw cap. Please indicate desired wire length (150 / 240 / 330 / 430 / 600 mm).
Measuring range from -200°C to +400°C.


Article No.	Probe Type	Thread Head	
40.440.10	1 x Pt-100 2-wire	SVL 15	
40.440.25	1 x Pt-100 4-wire	GL 18	
40.440.40	2 x Pt-100 2-wire	GL 18	

Portable Electronic Temperature Measuring Device TEST0 110

One channel meter for demanding temperature measurement in the laboratory, suitable for SCHMIZO 1 x Pt-100 4-wire probes

Article No.	Measuring Range	
40.441.00	- 100 ...+ 800 °C	

Battery and calibration protocol are included.
Input: USB Type C (use cable with LEMO - USB Type C plugs).

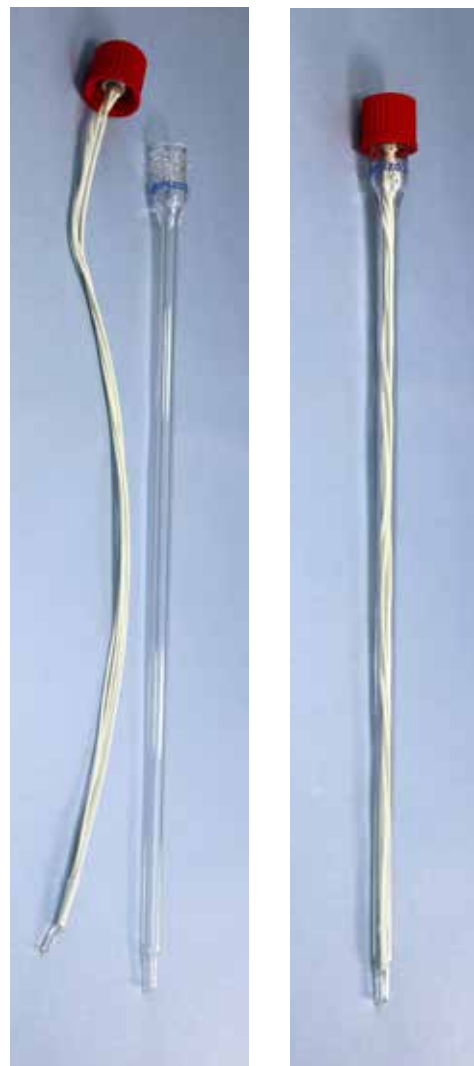

Probe Cables


Probe Cable for Pt-100

Cable length 200 cm, other length on request

Article No.	Plug	for Probe Type	Poles	
40.450.01	LEMO - none	1 x Pt-100 2-wire	2	
40.450.05	LEMO - none	1 x Pt-100 4-wire	4	
40.450.10	LEMO - none	2 x Pt-100 2-wire	4	
40.460.01	LEMO - LEMO	1 x Pt-100 2-wire	2	
40.460.05	LEMO - LEMO	1 x Pt-100 4-wire	4	
40.460.10	LEMO - LEMO	2 x Pt-100 2-wire	4	
40.465.50	LEMO - Mini-DIN	1 x Pt-100 4-wire	4	Testo 720
40.465.60	LEMO - USB-C	1 x Pt-100 4-wire	4	Testo 110



Probe Cable LEMO - LEMO

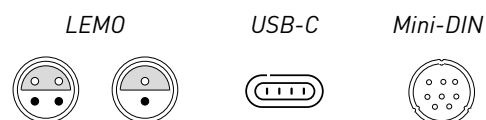


Probe Adapters

Guiding Tubes Ø 8 mm

for thermometers and probes; screw cap and joint are included

Art.Nr.	Ø mm	Form	NS	GL	
23.682.12	8	straight	14	18	
23.682.22	8	straight	19	18	
23.682.32	8	straight	24	18	
23.682.40	8	straight	29	18	
23.683.32	8	oblique	24	18	
23.683.40	8	oblique	29	18	

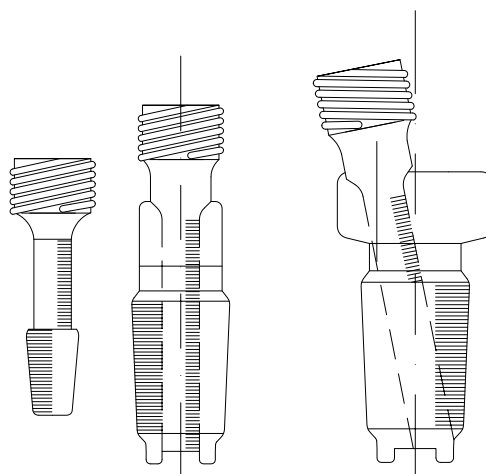

Probe adapter
straight

Probe adapter
EXCENTER
oblique

Spare Parts

for temperature probe adapters

Art.Nr.	Artikel	
12.255.18	Screw cap GL 18	
12.255.25	Screw cap GL 25	
12.263.08	Joint 16 x 8 mm	
12.265.12	Joint 22 x 12 mm	



pH and Redox Electrodes

Schmizo pH and Redox Electrodes

Low cost pH and redox electrodes without temperature compensation, for the full pH range from 0 - 14, applicable at temperatures from 0°C to 80°C.

Connection: SC plug head.

Together with the Schmizo guiding tube optimized for the use in our reactors.

Please specify desired usable length when ordering!

Article No.	Model	Total length mm	Usable length mm	
40.442.10	pH-Electrode	400 - 700	300 - 600	
40.442.11	Redox-Electr.	400 - 700	300 - 600	



Probe cable Schmizo pH-/Redox-Electrode

Ready made cable SC - BNC. Standard length 2 m (other lengths on request).

Article No.	Model	Ports	Length	
40.442.15	pH / Redox	SC - BNC	2 m	



Adapter for Schmizo pH-/Redox Electrodes

Guiding tube 12 mm straight, incl. cap GL 25 and seal 22x12 mm for pH-probes with shaft Ø12mm, NS cone 29/32

Article No.	Ø mm	Form	NS	GL	
23.682.45	12	straight	29	25	



Laboratory Racks

Table Racks for Reactors up to 5 Liters, Oblique Foot

made of anodized aluminium, tube Ø 34 mm made of stainless steel
2 cross connectors are included; height-adjustable feet

Article No.	Fits Reactor Holder		Dimensions w x h x d mm	
	Material	DN		
50.000.11	Polyester	60 / 100 / 120	500 x 500 x 1300	
50.000.16*	Polyester	150	500 x 500 x 1300	
50.000.21*	Stainl. Steel	120 / 150	500 x 500 x 1300	

* incl. 2 mounting rods

Art. 50.000.11 corresponds the picture on page 4

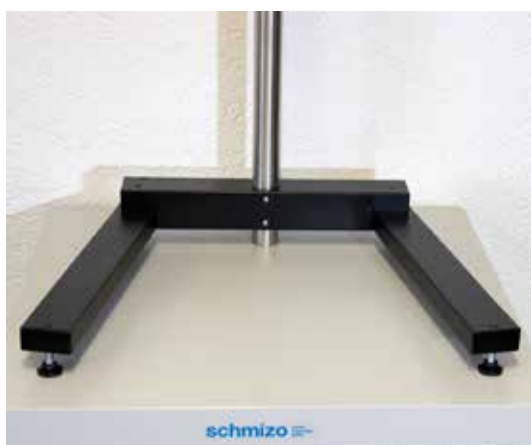


Table Racks for Reactors up to 5 Liters, Straight Foot


completely made of stainless steel, tube Ø 34 mm
2 cross connectors are included; height-adjustable feet

Article No.	Fits Reactor Holder		Dimensions w x h x d mm	
	Material	DN		
50.075.10	Polyester	60 / 100 / 120	350 x 500 x 1300	
50.075.45*	Polyester	150	350 x 500 x 1300	
50.075.55*	Stainl. Steel	120 / 150	350 x 500 x 1300	

* incl. 2 mounting rods



Oblique foot made of anodized aluminium

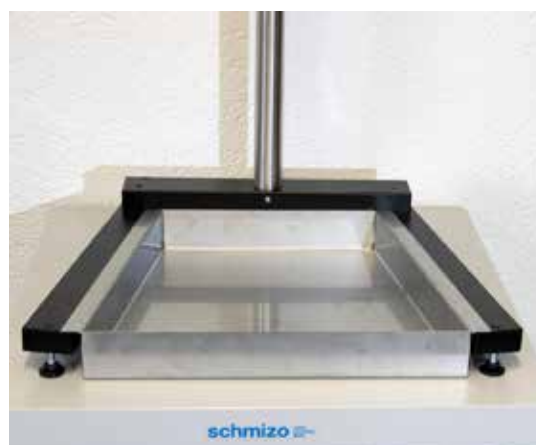

Oblique foot made with collecting tray made of stainless steel

Floor Standing Rack for Reactors up to 20 Liters

fitting reactor holders made of polyester DN 150 and reactor holders made of stainless steel DN 120 - 150; 1 cross connector is included. Strivings made of anodized aluminium, tube 1" (34 mm) made of stainless steel, height-adjustable feet or plastic rollers.

Article No.	Execution of the stand system	Dimensions w x h x d mm	
50.100.20	Adjustable feet	450 x 500 x 1300	
50.100.21	Plastic rollers	450 x 500 x 1500	



Floor Standing Rack for Reactors up to 20 Liters

fitting reactor holders DN 150 (polyester and stainless steel) and DN 200 (stainless steel); 1 cross connector is included. Completely made of stainless steel, tube 1" (34 mm), height-adjustable feet.

Article No.	Execution of the stand system	Dimensions w x d x h mm	
50.090.10	DN 150 / Plastic rollers	700 x 600 x 1800	
50.090.11	DN 150 / Adjustable feet	700 x 600 x 1800	
50.090.15	DN 200 / Plastic rollers	700 x 600 x 1800	
50.090.16	DN 200 / Adjustable feet	700 x 600 x 1800	



Art. 50.090.15 corresponds the picture on page 6

other sizes by request

Floor Standing Rack for Reactors up to 20 Liters with two Stackers

fitting reactor holders DN 150 (polyester and stainless steel) and DN 200 (stainless steel); 1 cross connector is included. Completely made of stainless steel, tube 1" (34 mm), height-adjustable feet or plastic rollers.

Article No.	Execution of the flange / stand system	Dimensions w x h x d mm	
50.080.10	DN 150 / Plastic rollers	700 x 690 x 1800	
50.080.11	DN 150 / Adjustable feet	700 x 690 x 1800	
50.080.15	DN 200 / Plastic rollers	700 x 690 x 1800	
50.080.16	DN 200 / Adjustable feet	700 x 690 x 1800	



All laboratory racks are also available with cross connectors fitting the **stirrer motors of the CAT 100-Series** (with 2 holders).


Cross Connector for CAT 100 Series Stirring Motors for 1 inch tubes

Article No.	Fixation	Ø1 mm	Ø2 mm	
50.010.40	Clamping screw	34	10 - 16	



Cross Connector Made of Anodized Aluminium for 1 inch tubes

Article No.	Fixation	Ø1 mm	Ø2 mm	
50.010.05	Allen Srew	34	10 - 16	
50.010.15	Clamping screw	34	10 - 16	

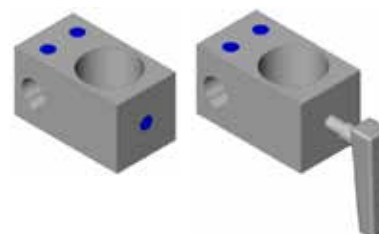

Cross Connector Made of Stainless Steel for 1 inch tubes

Article No.	Fixation	Ø1 mm	Ø2 mm	
50.120.15	Allen Srew	34	10 - 16	
50.120.20	Clamping screw	34	10 - 16	



Collecting Tray for Laboratory Rack made of stainless steel

Article No.	fits rack foot:	Volume l	Dimensions mm	
50.050.00	straight	6.5	400 x 280 x 70	
50.050.10	oblique	5.5	380 x 340 x 60	



Ground Rods for reactor holders

Article No.	Ø mm		
50.020.20	13	300	
50.020.25	13	490	
50.020.30	16	300	
50.020.35	16	490	
50.020.40	16	500	
50.020.45	16	550	

Supporting Ring silicone covered

Article No.	Ø mm	
50.050.50	13	
50.050.55	16	



Schlenk Glassware

We offer Schlenk lines in all possible variations. Since these lines are usually very individually assembled, we do not have fixed article numbers that we list here. We will be happy to make you an offer according to your needs.



Valves according to your choice:

- Stopcock with glass key
- Stopcock with PTFE key
- PRODURAN high vacuum valve
- YOUNG high vacuum valve




Schlenk line with YOUNG high vacuum valves

Schlenk Flasks

with NS socket and NS single or high vacuum valve

Article No.	Content	Neck NS	bore	
23.20x.39 *	25 ml	14/23	2.5	
23.20x.45 *	50 ml	14/23	2.5	
23.20x.55 *	100 ml	14/23	2.5	
23.20x.65 *	250 ml	14/23	2.5	
23.20x.55 *	100 ml	29/32	2.5	
23.20x.65 *	250 ml	29/32	2.5	
23.20x.73 *	500 ml	29/32	2.5	
23.20x.81 *	1000 ml	29/32	2.5	

* x = 0: valve with glass key, * x = 1: valve with PTFE key,
* x = 2: with PRODURAN valve, * x = 3: with YOUNG valve

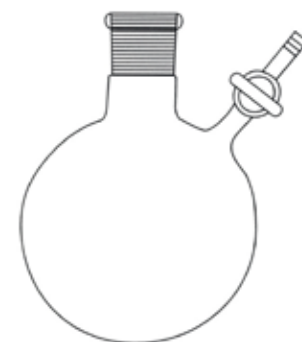


Schlenk Tubes

with NS socket and NS single or high vacuum valve

Article No.	Content	Neck NS	bore	
23.21x.31 *	10 ml	14/23	2.5	
23.21x.39 *	25 ml	14/23	2.5	
23.21x.45 *	50 ml	14/23	2.5	
23.21x.55 *	100 ml	14/23	2.5	
23.21x.65 *	250 ml	14/23	2.5	

All other sizes and combinations on demand

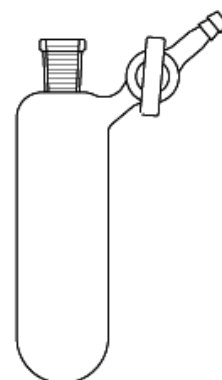


Products & Services Overview




Our reactor catalogue features the SCHMIZO reactors manufactured in our workshop and, of course, everything else you need to assemble a complete reactor system.



Our laboratory glassware catalogue features the current standard products and equipment (trade goods and products made in-house). And, if you don't see what you need here, we can get it for you.

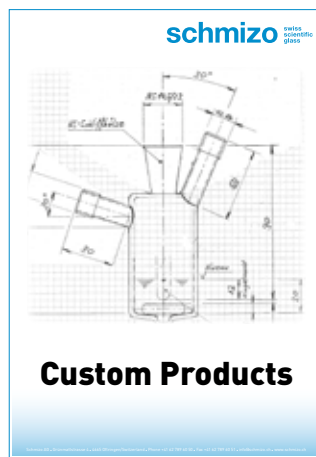


We specialise in manufacturing custom products according to your specifications, even as one-offs. Our technical consultant will gladly assist you in planning your project.

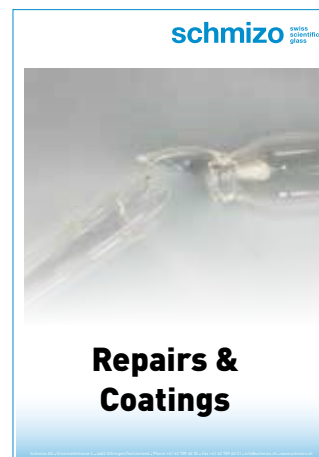

We perform glass repairs when they seem both reasonable and feasible. We also offer various coatings that protect against shattering and leaking in case of broken glassware.

Our Main Trading Partners

You can order these companies' products through Schmizo. We can also offer other brands if you wish.

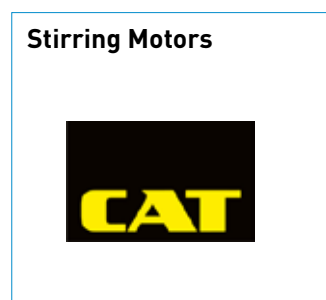

Technical Information

DURAN® Borosilicate Glass 3.3

For the manufacturing of our products we usually use DURAN® borosilicate glass 3.3. (manufactured by SCHOTT Rohrglas). The properties of DURAN® can be found on the homepage of DWK LIFE SCIENCES (www.dwk.com) or that of SCHOTT (www.schott.com) under the search term DURAN.

Pressure / Vacuum / Temperature

Schmizo reactors are designed for unpressurized operation. The authorized operating overpressure is generally 0.1 bar (but is dependent on the volume and type of vessel). Processed glass components in all nominal widths can be used at maximum vacuum, only vessels with flat bottom are not suitable for operation under slight overpressure and vacuum.

The maximum short-term operating temperature of borosilicate glass is 500° C, the maximum temperature shock must not exceed 100 Kelvin. Reactors with bottom valves made of PTFE can be operated at maximum temperature of 210° C. For jacketed reactors the maximum temperature difference ΔT between the reaction chamber and the heating jacket is 40 Kelvin and the maximum overpressure inside the heating jacket is 0.5 bar.

Before any evacuation or pressure load, glass vessels must be subjected a visual inspection for perfect condition. Damaged glass vessels should not be used to work with pressure and vacuum for safety reasons.

Coating of Glass Vessels

As splinter and leakage protection in case of glass breakage there are different kinds of coatings available. Gladly we make you a corresponding offer.

Special Designs

We are specialized in manufacturing customized system parts and can support you in development and realization with our specialist knowledge.

Glass Repairs

We perform all kinds of repairs if the defective component is made of borosilicate glass or quartz glass. In this case the breakage risk always remains by the contracting entity. For that we usually do not charge our costs in case of breakage.

Customer Service / Consultancy

We gladly stand at your disposal for a free consultation. Contact us by e-mail to info@schmizo.ch or by phone call to +41 (0)62 789 60 50.

Contacts

He looks for challenges

Ernst Schneeberger

Owner and CEO

ernst.schneeberger@schmizo.ch

+41(0)62 789 60 52



He always finds a solution

Sacha Ruch

Sales manager and Technical consultant

sacha.ruch@schmizo.ch

+41(0)62 789 60 55



He knows every product

Markus Lippuner

Sales and Purchases

markus.lippuner@schmizo.ch

+41(0)62 789 60 54



He makes your order running

Silas Huber

Sales Inside

silas.huber@schmizo.ch

+41(0)62 789 60 56



She knows how it works

Lisa De Lucia

Administration and Accounting

lisa.delucia@schmizo.ch

+41(0)62 789 60 53

